hyperSpec Plotting functions

Claudia Beleites <cbeleites@units.it> CENMAT, DMRN, University of Trieste

April 12, 2010

Reproducing the Examples in this Vignette

All spectra used in this manual are installed automatically with *hyperSpec*. Note that some definitions are executed in vignette.defs.

Contents

1	Predefined functions	1
2	Arguments for plot	3
3	Spectra 3.1 Stacked spectra	5 7
4	Levelplot	9
5	Spectra Matrix	9
6	Calibration Plots, (Depth) Profiles, and Time Series Plots	11
7	False-Colour Maps	11
8	3 D	13
9	Troubleshooting 9.1 No output is produced	13
10	Interactive Graphics 10.1 spc.identify: finding out wavelength, intensity and spectrum	13 13 13
	10.3 Related functions provided by base graphics and lattice	$\begin{array}{c} 14 \\ 14 \end{array}$

1 Predefined functions

hyperSpec comes with 5 major predefined plotting functions.

plot main switchyard for most plotting tasks


level plot hyperSpec has a level plot method

plotspc plots spectra

plotc calibration plot, time series, depth profile

plotmap version of levelplot allowing some more preprocessing of the data.


plotspc

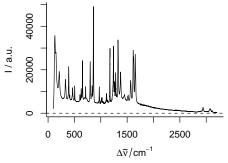

plots the spectra, i.e. the intensities \$spc over the wavelengths @wavelength.

> plotspc (paracetamol)

\mathbf{plotc}

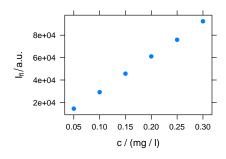


plots an intensity over a single other data column. e.g.

- time series
- calibration
- depth profile

> plotc (flu)

levelplot

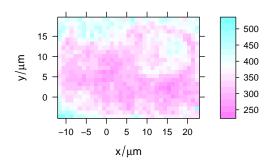


plots a false colour map, defined by a formula.

> levelplot (spc ~ x * y, chondro, aspect = "iso")

Warning: only the first wavelength is used.

plotmap

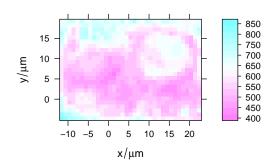

plotmap is a specialized version of level plot. It uses a single value (e.g. average intensity or cluster membership) over two data columns (default x and y

> plotmap (chondro)


2 Arguments for plot

The three specialized functions are also accessible via *plot*:

hyperSpec's plot method uses the second argument to determine which of the three specialized plot functions to call. All further arguments are handed over to this function.

This allows a few more handy abbreviations.


plot (x, "spc")

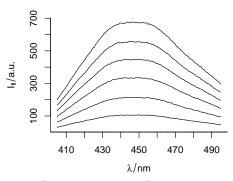

is equivalent to plotspc (flu)

> plot (flu, "spc")


plot (x, "spcmeansd")

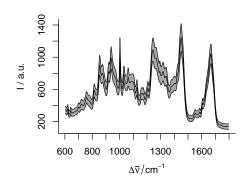


plots mean spectrum \pm 1 standard deviation

> plot (chondro, "spcmeansd")


plot (x, "spcprctile")

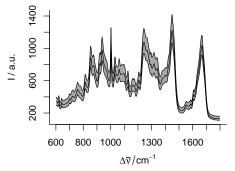


plot (x, "spcprctl5")

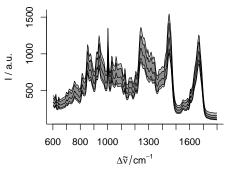

plot (x, "map")

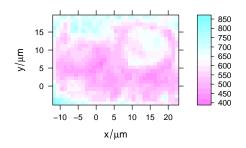

plot (x, "c")

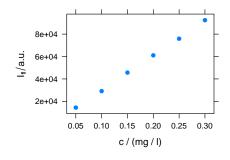

plots median, $16^{\rm th}$ and $84^{\rm th}$ percentile for each wavelength. For Gaussian distributed data, $16^{\rm th}$, $50^{\rm th}$ and $84^{\rm th}$ percentile are equal to mean \pm standard deviation. Spectroscopic data frequently are not Gaussian distributed. The percentiles give a better idea of the true distribution. They are also less sensitive to outliers.

> plot (chondro, "spcprctile")

like "spcprct1" plus $5^{\rm th}$ and $95^{\rm th}$ percentile.

> plot (chondro, "spcprct15")


is equivalent to plotmap (chondro)

> plot (chondro, "map")

is equivalent to plotc (flu)

> plot (flu, "c")


plot (x, "ts")

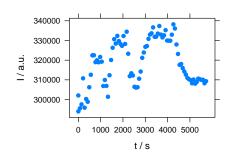

plots a time series plot, equivalent to plotc (laser, spc $\tilde{\ }$ t)

> plot (laser, "ts")

plot (x, "depth")

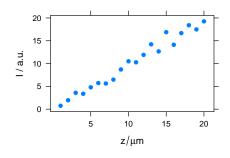


plots a depth profile plot, equivalent to plotc (laser, spc \tilde{z})

```
> depth.profile <- new ("hyperSpec",
+     spc = as.matrix (rnorm (20) + 1:20),
+     data = data.frame (z = 1 : 20),
+     label = list (spc = "I / a.u.",
+     z = expression (`/` (z, mu*m)),
+     .wavelength = expression (lambda)))</pre>
```

> plot (depth.profile, "depth")

plot (x, "mat")

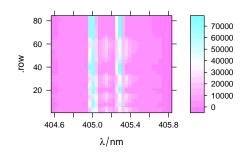

plots the spectra matrix.

> plot (laser, "mat")

Equivalent to


> levelplot (spc ~ .wavelength * .row, laser)

3 Spectra


plotspc

plotspc offers a variety of parameters for customized plots. To plot ...


particular wavelength ranges

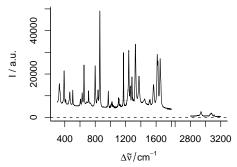


particular wavelength ranges II

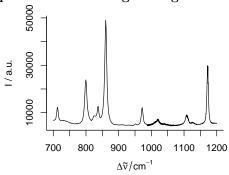

with reversed abscissa

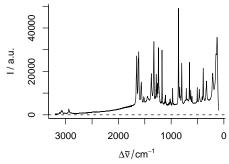

in different colours

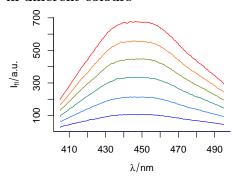

use wl.range = list (600 ~ 1800, 2800 ~ 3100). If wl.range already contains indices: use wl.index = TRUE. Cut the wavelength axis appropriately with xoffset = 750

```
> plotspc (paracetamol,
+ wl.range = c (300 ~ 1800, 2800 ~ max),
+ xoffset = 750)
```


if only one wavelength range is needed, the extract command is handier:

> plotspc (paracetamol[,, 700 ~ 1200])

```
use wl.reverse = TRUE
> plotspc (paracetamol, wl.reverse = TRUE )
```


```
use col = vector.of.colours
> plotspc (flu, col = matlab.dark.palette(6) )
```

dots instead of lines

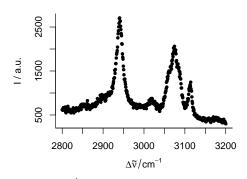


```
use lines.args = list (pch = 20, type = "p")
> plotspc (paracetamol [,, 2800 ~ 3200],
+ lines.args = list (pch = 20, type = "p"))
```

${f mass}$ spectra

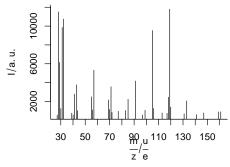

use lines.args = list (type = "h")
> plot (barbituates [[1]], lines.args = list (type = "h"))


more spectra into an exsisting plot

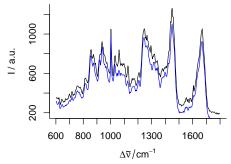

use add = TRUE

> plotspc (chondro[30,,])
> plotspc (chondro[300,,], add = TRUE, col = "blue")

with different line at I=0

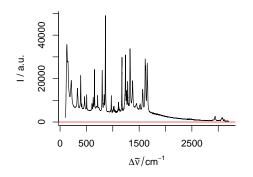


 $\begin{tabular}{ll} use & {\tt zeroline = list.of.arguments.to.abline.} \\ {\tt NULL suppresses the line.} \\ \end{tabular}$

> plotspc (paracetamol,
+ zeroline = list (col = "red"))

3.1 Stacked spectra


stacked = TRUE

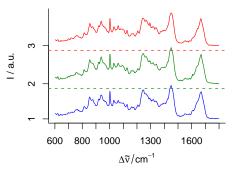

use stacked = TRUE


```
> plotspc (chondro[1:3,,],
+ col = matlab.dark.palette (3),
+ stacked = TRUE)
```

Stacking groups of spectra

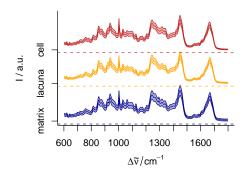


if you want to plot mean spectrum \pm 1 standard deviation of each of the clusters, use stacked = grouping.column ".name"

yoffset

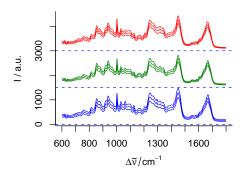

or use yoffset =

Manually with normalized spectra

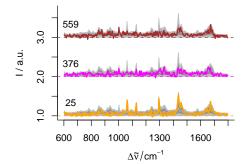

```
it's also possible to normalize and stack the spectra like this:  \\
```

```
> out <- c (25, 376, 559)
> outcols <- c ("orange", "magenta", "brown")</pre>
Preprocessing of the spectra:
> bl <- spc.fit.poly.below (chondro)
Fitting with npts.min = 15
> chondro <- chondro - bl
> chondro <- sweep (chondro, 1, mean, "/")</pre>
 chondro <- sweep (chondro,
                     2,
                     apply (chondro,
                            2,
                            quantile,
                            0.05),
The figure:
> ## coordinate system:
> plot(chondro[1],
       plot.args = list (ylim = c (1, length (out) + .7)),
       lines.args = list( type = "n")
> ## stacked spectra:
 for (i in seq (along = out)){
    plot(chondro,
         "spcprct15",
         yoffset = i,
         col = "gray",
         add = TRUE)
    plot (chondro [out[i]],
          yoffset = i,
          col = outcols[i],
          add = TRUE,
          lines.args = list (lwd = 2))
```

text (650, i + .33, out [i]) }

4 Levelplot

Levelplot can use two special column names:

```
.wavelength for the wavelengths
```


.row for the row index (i.e. spectrum number) in the data

5 Spectra Matrix


It is often useful to plot the spectra against an additional coordinate, e.g. the time for time series, the depth for depth profiles, etc.

This can be done by plot (object, "mat") or levelplot (model = spc ~.wavelength * other.data.column, object). The actual plotting is done by levelplot, so the plots can be grouped or conditioned.

different palette

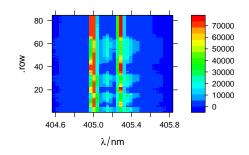


different y axis

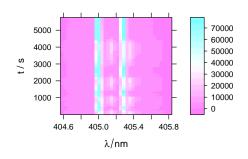

Changing the y axis is only possible with levelplot:

```
> levelplot (spc ~ .wavelength * t,
+ laser)
```


contour lines with plot

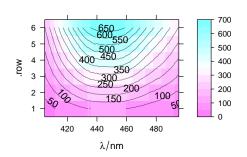

Contourplots are possible with plot and levelplot:

```
> plot (flu,
+          "mat",
+          contour = TRUE,
+          labels = TRUE,
+          col = "#00000080",
+          at = seq (0, 700, by = 50))
```

contour lines with levelplot

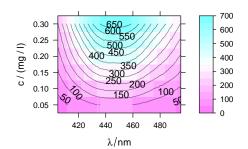


Contourplot with a different y axis:


```
> levelplot (spc ~ .wavelength * c,
+ flu,
+ contour = TRUE,
+ labels = TRUE,
+ col = "#00000080",
+ at = seq (0, 700, by = 50))
```

6 Calibration Plots, (Depth) Profiles, and Time Series Plots

plotc

plotc plots an intensity over one of the extra data columns. The abscissa uses column \$c\$ by default, another column can be specified using a proper formula. The ordinate can be calculated as a sum characteristic (with parameter func= function, defaulting to sum).

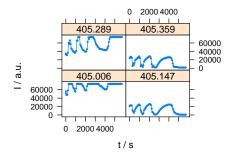

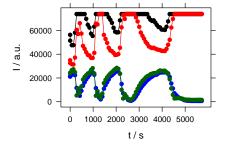

7 False-Colour Maps


plotmap

plotmap uses levelplot, a *lattice* function. Therefore, in loops, functions, Sweave chunks, etc. the lattice object needs to be printed explicitly by print (plotmap (object)) (R FAQ: Why do lattice/trellis graphics not work?).

plotmap produces a 3d plot, with the z axis colour-coded. plotmap's arguments x and y take the name of extra data columns.


The colour-coded axis. Also z can be used to select one column of the extra data by name. Alternatively, it may be a numeric or factor directly giving the values to be used. Each level of a factor will have one colour. It is also possible to plot a sum characteristic of the spectra: supply the function in argument func. The default setting is to plot the average intensity (no z and func=mean).

plotting clusters

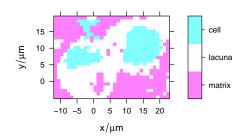


> plotmap (chondro, clusters ~ x * y)

different palette

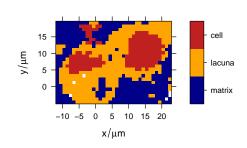

To plot with a different palette, use *trellis.args=* list (col.regions = palette).

defined wavelengths

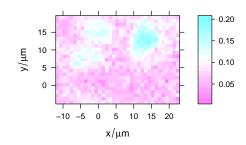

To plot a map with particular wavelengths use this:

```
> plotmap (chondro[, , c( 728, 782, 1098, 1240,
+ 1482, 1577)])
```

Examples

Conditioning. Lattice graphics have a concept of conditioning a plot. Instead of plotting all data in one diagram, a diagram is produced for each of the groups specified by the condition. plotmap's argument *cond* takes he name of the extra data column used for conditioning. This could e.g. be a column containing the sample number of a *hyperSpec* object that contains several samples.

Beispiele: voronoi

8 3 D


9 Troubleshooting

9.1 No output is produced

plotmap and plotc use levelplot, a *lattice* function. Therefore, in loops, functions, Sweave chunks, etc. the lattice object needs to be printed explicitly by print (plotmap (object)) (R FAQ: Why do lattice/trellis graphics not work?).

In vignettes.defs is shown how the lattice functions can be redefined so that the result is printed without external print command.

in functions

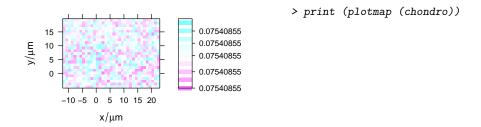

10 Interactive Graphics

hyperSpec offers two basic interaction functions, spc.identify, and map.identify. They identify points in spectra plots and map plots, respectively.

$10.1~{ m spc.identify:}$ finding out wavelength, intensity and spectrum

spc.identify allows to measure points in graphics produced by plotspc. It works correctly with reversed and cut wavelength axes.

```
> spc.identify (plotspc (paracetamol, wl.range = c (600 ~ 1800, 2800 ~ 3200), xoffset = 800))
```

The result is a data.frame with the indices of the spectra, the wavelength, and its intensity.

10.2 map.identify: finding a spectrum in a map plot

map.identify returns the spectra indices of the clicked points.

> map.identify (chondro)

10.3 Related functions provided by base graphics and lattice

For base graphics (as produced by plotspc), locator may be useful as well. It returns the clicked coordinates. Note that these are *not* transformed according to xoffset & Co.

For lattice graphics, grid.locator may be used instead. If it is not called in the panel function, a preceeding call to trellis.focus is needed:

```
> plot (laser, "mat")
> trellis.focus ()
> grid.locator ()
```

identify (or panel.identify for lattice graphics) allows to identify points of the plot directly. Note that the returned indices correspond to the plotted object.

10.4 Interactively changing graphics

hyperSpec's lattice functions work with playwith and latticist. These packages allow easy customization of the plots and also identification of points.