

# Calibration of Quinine Fluorescence Emission Vignette for the Data Set `flu` of the R package `hyperSpec`

Claudia Beleites ([cbeleites@units.it](mailto:cbeleites@units.it))  
CENMAT, DMRN, University of Trieste

December 14, 2009

## Reproducing the Examples in this Vignette

The spectra files are shipped with *hyperSpec*. This allows reproduction of the whole vignette (source code and spectra files are in the documentation directory).

For reproducing the examples in a live session, the full file names of the spectra can be found with the command:

```
> Sys.glob (paste (.libPaths (), "hyperSpec/doc/flu*.txt", sep = "/"))
```

This vignette gives an example how to

- write an import function for a spectrometer manufacturer's proprietary ASCII files,
- add further data columns to the spectra, and
- set up a linear calibration (inverse least squares).

The data set `flu` in *hyperSpec* consists of 6 fluorescence emission spectra of quinine solutions. They were acquired during an student practicum and were kindly provided by M. Kammer.

The concentrations of the solutions range from 0.05 mg/l to 0.30 mg/l. Spectra were acquired with a Perkin Elmer LS50-B fluorescence spectrometer at 350 nm excitation.

In order to work with *hyperSpec*, the package needs to be loaded in R:

```
> library (graphics)  
> library (hyperSpec)
```

## 1 Writing an Import Function

Now we need to import the spectra. They are in Perkin Elmer's ASCII file format, one spectrum per file. The files are completely ASCII text, with the actual spectra starting at line 55.

The function should automatically read in all files specified by a pattern, such as `*.txt`. In order to gain some speed, the spectra matrix is preallocated after the first file is read. Also, the number of header lines are predefined instead of searching for the line after `#DATA`.

Note, that labels giving the correct units (e.g. for axis labels) are set. The label with the special name `.wavelength` corresponds to the wavelength axis, all data columns should have a label with the same name. The spectra are always in a data column called `spc`.

```
> read.PE <- function (files = "flu?.txt", skip = 54) {
+   files <- Sys.glob (files)
+
+   buffer <- matrix (scan (files [1], skip = skip), ncol = 2, byrow = TRUE)
+
+   wavelength <- buffer [, 1]
+   spc <- matrix (ncol = nrow (buffer), nrow = length (files))
+
+   spc [1, ] <- buffer [, 2]
+
+   for (f in seq (along = files)[-1]) {
+
+       buffer <- matrix (scan (files [f], skip = skip), ncol = 2, byrow = TRUE)
+
+       if (! all.equal (buffer [, 1], wavelength))
+           stop (paste(files [f], "has different wavelength axis."))
+
+       spc [f, ] <- buffer[, 2]
+
+   }
+
+   new ("hyperSpec", wavelength = wavelength, spc = spc,
+       label = list (.wavelength = expression (lambda[f1] / nm),
+           spc = "I / a.u."))
+ }
```

From now on, the function can be used:

```
> flu <- read.PE ("flu?.txt")
```

Now the spectra are in a *hyperSpec* object and can be examined e.g. by

```
> flu
```

hyperSpec object

6 spectra

1 data columns


181 data points / spectrum

wavelength: lambda[f1]/nm [numeric 181] 405.0 405.5 ... 495

data: (6 rows x 1 columns)

(1) spc: I / a.u. [AsIs matrix 6 x 181] range 27.15000 32.34467 ... 677.4947

```
> plot (flu)
```

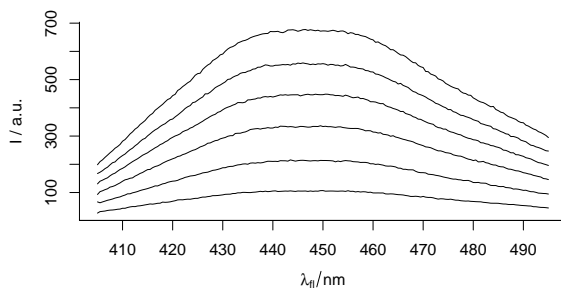

## 2 Adding further Data Columns

The calibration model needs the quinine concentrations for the spectra. This information can be stored together with the spectra, and also gets an appropriate label:

```
> flu$c <- seq (from = 0.05, to = 0.30, by = 0.05)
> labels (flu, "c") <- "c / (mg / l)"
> flu

hyperSpec object
  6 spectra
  2 data columns
  181 data points / spectrum
wavelength: lambda[f1]/nm [numeric 181] 405.0 405.5 ... 495
data: (6 rows x 2 columns)
  (1) spc: I / a.u. [AsIs matrix 6 x 181] range 27.15000 32.34467 ... 677.4947
  (2) c: c / (mg / l) [numeric 6] range 0.05 0.10 ... 0.3
```

Now the *hyperSpec* object *flu* contains two data columns, holding the actual spectra and the respective concentrations. The dollar operator returns such a data column:

```
> flu$c


[1] 0.05 0.10 0.15 0.20 0.25 0.30
```

## 3 Linear Calibration

As R is developed for the purpose of statistical analysis, tools for a least squares calibration model are readily available.

The original spectra range from 405 to 495 nm. However, the intensities at 445 nm are perfect for a univariate calibration. Plotting them over the concentration is done by:

```
> plotc (flu[, ,445])
```

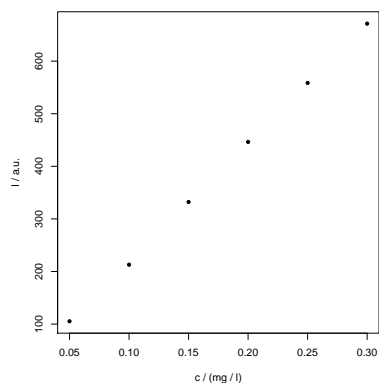


The square bracket operator extracts parts of a *hyperSpec* object. The first coordinate defines which spectra are to be used, the second which data columns, and the third gives the spectral range.

We discard all the wavelengths but 445 nm:

```
> flu <- flu[,445]
```

The plot could be enhanced by annotating the ordinate with the emission wavelength. Also the axes should start at the origin, so that it is easier to see whether the calibration function will go through the origin:

```
> plotc(flu, zlab = expression(I["450 nm"] / a.u.),
+       plot.args = list(xlim = range(0, flu$c), ylim = range(0, flu$spc)))
```

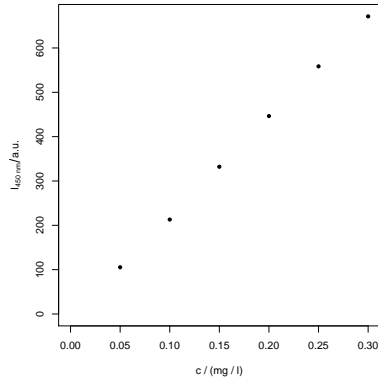

The actual calibration is a linear model, which can be fitted by the R function `lm`. `lm` needs a *formula* that specifies which data columns are dependent and independent variables.

The normal calibration plot gives the emission intensity as a function of the concentration, and the calibration function thus models  $I = f(c)$ , i. e.  $I = mc + b$  for a linear calibration. This is then solved for  $c$  when the calibration is used.

However, R's linear model is a quite strict in predicting: a model set up as  $I = f(c)$  will predict the intensity as a function of the concentration but not the other way round. Thus we set up an inverse calibration model<sup>1</sup>:  $c = f(I)$ . The corresponding R formula is  $c \sim I$ , or in our case  $c \sim \text{spc}$ , as the intensities are stored in the data column `$spc`:

In addition, `lm` (like most R model building functions) expects the data to be a *data.frame*.

There are three abbreviations that help to get the parts of the *hyperSpec* object that are frequently needed:

<code>flu[[ ]]</code>	returns the spectra matrix. It takes the same indices as <code>[ ]</code> .
<code>flu\$.</code>	returns the data as a <i>data.frame</i>
<code>flu\$..</code>	returns a <i>data.frame</i> that has all data columns but the spectra

```
> flu[[ ]]
```

```
      445
```

```
[1,] 105.3617
```

```
[2,] 212.9753
```

```
[3,] 332.1467
```

```
[4,] 446.3653
```

```
[5,] 558.4950
```

```
[6,] 671.2343
```

---

<sup>1</sup>As we can safely assume that the error on the concentrations is far larger than the error on the instrument signal, it is actually the correct type of model from the least squares fit point of view.

```
> flu$.
      445      c
1 105.361666 0.05
2 212.975333 0.10
3 332.146666 0.15
4 446.365333 0.20
5   558.495 0.25
6 671.234333 0.30
```

```
> flu$..
```

```
      c
1 0.05
2 0.10
3 0.15
4 0.20
5 0.25
6 0.30
```

Putting this together, the calibration model is calculated:

```
> calibration <- lm (c ~ spc, data = flu$.)
```

The `summary` gives a good overview of our model:

```
> summary (calibration)
```

Call:

```
lm(formula = c ~ spc, data = flu$.)
```

Residuals:

```
      1      2      3      4      5      6
-0.0008475  0.0018421 -0.0005493 -0.0007634 -0.0000591  0.0003772
```

Coefficients:

```
      Estimate Std. Error t value Pr(>|t|)
(Intercept) 4.527e-03  1.034e-03   4.38  0.0119 *
spc         4.396e-04  2.383e-06 184.47 5.18e-09 ***
---

```

```
Signif. codes:  0 '***' 0.001 '**' 0.01 '*' 0.05 '.' 0.1 ' ' 1
```

Residual standard error: 0.001134 on 4 degrees of freedom

Multiple R-squared: 0.9999, Adjusted R-squared: 0.9999

F-statistic: 3.403e+04 on 1 and 4 DF, p-value: 5.18e-09

In order to get predictions for new measurements, a new *data.frame* with the same independent variables (in columns with the same names) as in the calibration data are needed. Then the function `predict` can be used. It can also calculate the prediction interval. If we observe e.g. an intensity of 125 units, the corresponding concentration and its 99 % prediction interval are:

```
> I <- 125
> conc <- predict (calibration, newdata = list (spc = as.matrix(I)), interval = "prediction",
+               level = .99)
> conc
```

	fit	lwr	upr
1	0.05948115	0.05314837	0.06581394

Finally, we can draw the calibration function and its 99 % confidence interval (also via `predict`) together with the prediction example:

```
> plotc (flu, zlab = expression (I ["450 nm"] / a.u.),
+       plot.args = list (xlim = range (0, flu$c), ylim = range (0, flu$spc)))
> int <- list (spc = as.matrix(min (flu) : max(flu)))
> ci <- predict (calibration, newdata = int, interval = "confidence", level = 0.99)
> matlines (ci, int$spc, col = c ("red", "gray", "gray"), lty = 1)
> # extrapolate to lower intensities
> int <- list (spc = as.matrix(0 : min (flu)))
> ci <- predict (calibration, newdata = int, interval = "confidence", level = 0.99)
> matlines (ci, int$spc, col = c ("red", "#606060", "#606060"), lty = 3)
> # our example
> lines (conc[-1], rep(I, 2), col = "blue")
> points (conc[1], I, col = "blue", pch = 4, cex = 0.5)
```

