

# hyperSpec Introduction

Claudia Beleites <[cbeleites@units.it](mailto:cbeleites@units.it)>  
CENMAT, DMRN, University of Trieste

February 13, 2011

## Reproducing the Examples in this Vignette

All spectra used in this manual are installed automatically with *hyperSpec*.  
Note that some definitions are executed in `vignette.defs`.

## Contents

<b>1. Introduction</b>	<b>3</b>
1.1. Notation . . . . .	3
<b>2. Remarks on R</b>	<b>4</b>
2.1. Generic Functions . . . . .	4
2.2. Functionality Can be Extended at Runtime . . . . .	4
2.3. Validity Checking . . . . .	4
2.4. Special Function Names . . . . .	4
2.4.1. The Names of Operators . . . . .	4
2.4.2. Assignment Functions . . . . .	5
<b>3. Loading and the package and configuration</b>	<b>5</b>
<b>4. The structure of hyperSpec objects</b>	<b>5</b>
<b>5. Functions provided by hyperSpec</b>	<b>6</b>
<b>6. Obtaining Basic Information about hyperSpec Objects</b>	<b>6</b>
<b>7. Creating a hyperSpec Object, Data Import and Export</b>	<b>7</b>
7.1. Creating a <i>hyperSpec</i> Object from Spectra Matrix and Wavelength Vector . . . . .	7
<b>8. The Logbook</b>	<b>7</b>
<b>9. Combining and Decomposing hyperspec Objects</b>	<b>8</b>
9.1. Binding Objects together . . . . .	8
9.2. Binding Objects that do not Share the Same Extra Data and/or Wavelength Axis . . . . .	9
9.3. Binding Objects that do not Share the Same Spectra . . . . .	9
9.4. Matrix Multiplication . . . . .	11
9.5. Decomposition . . . . .	11

<b>10. Access to the data</b>	<b>11</b>
10.1. Selecting and Deleting Spectra	11
10.1.1. Random Samples	12
10.1.2. Sequences	13
10.2. Selecting Extra Data Columns	13
10.3. Selecting Wavelength Ranges	17
10.4. Deleting Wavelength Ranges	18
10.4.1. Converting Wavelengths to Indices and vice versa	19
10.4.2. Changing the Wavelength Axis	20
10.4.3. Ordering the Wavelength Axis	21
10.5. More on the Square-Bracket Operators for Replacing Values	21
10.6. Fast Access to Parts of the <i>hyperSpec</i> Object	22
10.7. Conversion to Long-Format <i>data.frame</i>	22
<b>11. Plotting</b>	<b>22</b>
<b>12. Spectral (Pre)processing</b>	<b>22</b>
12.1. Cutting the Spectral Range	22
12.2. Shifting Spectra	23
12.2.1. Calculating the Shift	23
12.3. Smoothing Interpolation	24
12.4. Background Correction	26
12.5. Offset Correction	26
12.6. Baseline Correction	26
12.7. Intensity Calibration	26
12.7.1. Correcting by a constant, e. g. Readout Bias	26
12.7.2. Correcting Wavelength Dependence	26
12.7.3. Spectra Dependent Correction	27
12.8. Normalization	27
12.9. Centering the Data	27
12.10 Variance Scaling	28
12.11 Multiplicative Scatter Correction (MSC)	28
12.12 Spectral Arithmetic	28
<b>13. Data Analysis</b>	<b>29</b>
13.1. Data Analysis Methods using a <i>data.frame</i>	
e. g. Principal Component Analysis with <i>prcomp</i>	29
13.1.1. PCA as Noise Filter	30
13.2. Data Analysis using long-format <i>data.frame</i>	
e. g. plotting with <i>ggplot2</i>	31
13.3. Data Analysis Methods using a matrix	
e. g. Hierarchical Cluster Analysis	31
13.4. Calculating group-wise Sum Characteristics	
e. g. Cluster Mean Spectra	31
13.5. Splitting an Object, and Binding a List of <i>hyperSpec</i> Objects	31
<b>14. Speed Considerations</b>	<b>33</b>
<b>A. Overview of the functions provided by <i>hyperSpec</i></b>	<b>33</b>

## 1. Introduction

*hyperSpec* is a R package that allows convenient handling of (hyper)spectral data sets, i. e. data sets comprising spectra together with further data on a per-spectrum basis. The spectra can be anything that is recorded over a common discretized axis.

*hyperSpec* works with any data that fits in this general scheme, so that the three terms may also be used for:

wavelength: frequency, wavenumbers, chemical shift, Raman shift,  $\frac{m}{z}$ , etc.

intensity: transmission, absorbance,  $\frac{e^-}{s}$ , ...

extra data: spatial information (spectral images, maps, or profiles), temporal information (kinetics, time series), concentrations (calibration series), class membership information, etc.  
Note that there is no restriction on the number of extra data columns.

Throughout the documentation of the package, the terms intensity and wavelength axes refer to the spectral ordinate and abscissa, respectively.

This vignette gives an introduction on basic working techniques using the R package *hyperSpec*.

*hyperSpec* comes with five data sets,

**chondro** a Raman map of chondrocytes in cartilage,

**flu** a set of fluorescence spectra of a calibration series, and

**laser** a time series of an unstable laser emission

**paracetamol** a Raman spectrum of paracetamol (acetaminophene) ranging from 100 to 3200  $\text{cm}^{-1}$  with overlapping wavelength ranges.

**barbituates** GC-MS spectra with differing wavelength axes as a list of 286 *hyperSpec* objects.

In this vignette, the data sets are used to illustrate appropriate procedures for different tasks and different spectra. In addition, the first three data sets are accompanied by their own vignettes showing exemplary work flows for the respective data type.

This document describes how to accomplish typical tasks in the analysis of spectra. It does not give a complete reference on particular functions. It is therefore recommended to look up the methods in R's help system using `? command`.

A list of all functions available in *hyperSpec* is given in appendix A (p. 33).

### 1.1. Notation

This vignette demonstrates working techniques mostly from a spectroscopic point of view: rather than going through the functions provided by *hyperSpec*, it is organized more closely on spectroscopic tasks. However, the functions discussed are printed on the margin for a fast overview.

In R, slots of a S4 class can be accessed directly by the `@` operator. In this vignette, the notation `@xxx` will thus mean “slot *xxx* of an object”. Likewise, named elements of a *list*, like the columns of a *data.frame*, are accessed by the `$` operator, and `$xxx` will be used for “column *xxx*”, and as an abbreviation for “column *xxx* of the *data.frame* in slot *data* of the object”.

## 2. Remarks on R

### 2.1. Generic Functions

*Generic Functions* are functions that apply to a wide range of data types or classes, e.g. *plot*, *print*, mathematical operators, etc. These functions can be implemented in a specialized way by each class. *hyperSpec* implements with a variety of such functions, see the table in appendix A on page 33.

### 2.2. Functionality Can be Extended at Runtime

R's concept of functions offers much flexibility. Functions may be added or changed by the user in his *workspace* at any time. This is also true for methods belonging to a certain class. Neither restart of R nor reloading of the package or anything the like is needed. If the original function resides in a namespace (as it is the case for all functions in *hyperSpec*), the original function is not deleted. It is just masked by the user's new function but stays accessible via the `::` operator.

This offers the opportunity of easily writing specialized functions that are adapted to specific tasks. *hyperSpec*'s vignettes use this to set up special versions of the lattice graphics functions that are already wrapped in `print` (see also [R FAQ: Why do lattice/trellis graphics not work?](#)) and allow the code in the code chunks of the vignettes to be exactly what one would type during an interactive R session. For the code, check the `vignettes.defs` file accompanying all *hyperSpec* vignettes.

### 2.3. Validity Checking

S4 classes have a mechanism to define and enforce that the data actually stored in the object is appropriate for this class. In other words, there is a mechanism of *validity checking*.

The functions provided by *hyperSpec* check the validity of *hyperSpec* objects at the beginning, and – if the validity could be broken by inappropriate arguments – also before leaving the function.

It is highly recommended to use validity checking also for user-defined functions. In addition, non-generic functions should first ensure that the argument actually is a *hyperSpec* object. The two tasks are accomplished by:

```
> chk.hy (object)
> validObject (object)
```

The first line checks whether `object` is a *hyperSpec* object, the second checks its validity. Both functions return `TRUE` if the checks succeed, otherwise they raise an error and stop.

### 2.4. Special Function Names

#### 2.4.1. The Names of Operators

Operators such as `+`, `*`, `%`, etc. are in fact functions in R. Thus they can be handed over as arguments to other functions (particularly to the vectorization functions `*apply`, `sweep`, etc.). In this case the name of the function must be quoted: ``*`` is the recommended style (although `"*"` will often work as well), e.g.:

```
> sweep (flu, 2, mean, ``*)
```

These functions can also be called in a more function-like style (prefix notation):

```
> `+` (3, 5)
[1] 8
```

slot	get	set
@wavelength	wl	wl<-
@data	[, [[, \$, as.data.frame, as.long.df, ...	[<-, [[<-, \$<-
@label	labels	labels<-
@log	logbook	logentry

Table 1: Get and set functions for the slots of *hyperSpec* objects

### 2.4.2. Assignment Functions

R allows the definition of functions that do an assignment (set some part of the object), such as:

```
> wl (flu) <- new.wavelength.values
```

an assignment to variable `wl`: `wl<-``.

## 3. Loading and the package and configuration

To load *hyperSpec*, use

```
> library(hyperSpec)
```

The global behaviour of *hyperSpec* can be configured via options. The values of the options are retrieved with `hy.getOptions` and `hygetOption`, and changed with `hy.setOptions`.

Currently, the only option provided is `log`, a logical specifying whether assignment functions should automatically add entries to the logbook (see section 8, p. 7).

## 4. The structure of *hyperSpec* objects

*hyperSpec* is a S4 (or new-style) class. Four slots contain the parts of the object:

**@wavelength** containing a numeric vector with the wavelength axis of the spectra.

**@data** a *data.frame* with the spectra and all further information belonging to the spectra

**@label** a list with appropriate labels (particularly for axis annotations)

**@log** a *data.frame* keeping track of what is done with the object

While the parts of the *hyperSpec* object can be accessed directly, it is good practice to use the functions provided by *hyperSpec* to handle the objects rather than accessing the slots directly (tab. 1). This also ensures that proper (*valid*) objects are retained.

Most of the data is stored in **@data**. This *data.frame* has one special column, **\$spc**. It is the column that actually contains the spectra. The spectra are stored in a matrix inside this column, as illustrated in figure 1. Even if there are no spectra, **\$spc** must still be present. It is then a matrix with zero columns.

Slot **@label** contains an element for each of the columns in **@data** plus one holding the label for the wavelength axis, **.wavelength**. The elements of the list may be anything suitable for axis annotations, i.e. they should be either character strings or expressions for “pretty” axis annotations (see e.g. figure 5 on page 25). To get familiar with expressions for axis annotation, see `?plotmath` and `demo(plotmath)`.

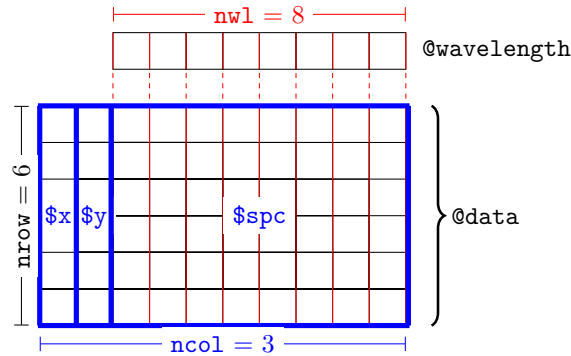


Figure 1: The structure of the data in a *hyperSpec* object.

## 5. Functions provided by hyperSpec

Table A (p. 33) in the appendix gives an overview of the functions implemented by *hyperSpec*.

## 6. Obtaining Basic Information about hyperSpec Objects

As usual, the *print* and *show* methods display information about the object, and *summary* yields some additional details about the data handling done so far:

*print*, *show*,  
*summary*

```
> chondro

hyperSpec object
  875 spectra
  4 data columns
  300 data points / spectrum
wavelength: Delta * tilde(nu)/cm^-1 [numeric] 602 606 ... 1798
data: (875 rows x 4 columns)
  1. y: y/(mu * m) [numeric] -4.77 -4.77 ... 19.23
  2. x: x/(mu * m) [numeric] -11.55 -10.55 ... 22.45
  3. clusters: clusters [factor] matrix matrix ... lacuna + NA
  4. spc: I / a.u. [matrix300] 501.82 500.46 ... 169.29

> summary (chondro)

hyperSpec object
  875 spectra
  4 data columns
  300 data points / spectrum
wavelength: Delta * tilde(nu)/cm^-1 [numeric] 602 606 ... 1798
data: (875 rows x 4 columns)
  1. y: y/(mu * m) [numeric] -4.77 -4.77 ... 19.23
  2. x: x/(mu * m) [numeric] -11.55 -10.55 ... 22.45
  3. clusters: clusters [factor] matrix matrix ... lacuna + NA
  4. spc: I / a.u. [matrix300] 501.82 500.46 ... 169.29
log:
      short      long      date      user
1 .local  examp...  2011-02-13 22:10:51  cb@cb
```

The data set *chondro* consists of 875 spectra with 300 data points each, and 4 data columns: two for the spatial information, one factor with the results of a cluster analysis plus *\$spc*. These informations can be directly obtained by

*nrow*, *ncol*,  
*nwl*, *dim*

```
> nrow (chondro)
```

```
[1] 875

> nwl (chondro)

[1] 300

> ncol (chondro)

[1] 4

> dim (chondro)

nrow ncol nwl
875    4  300
```

The names of the columns in `@data` are accessed by

```
> colnames (chondro)

[1] "y"          "x"          "clusters" "spc"
```

```
colnames,
rownames,
dimnames, w1
```

Likewise, `rownames` returns the names assigned to the spectra, and `dimnames` yields a list of these three vectors (including also the column names of `$spc`). The column names of the spectra matrix are the wavelengths. They are accessed by `w1`, see section [10.4.2](#).

Extra data column names and rownames of the object may be set by `colnames<-` and `rownames<-`, respectively. `colnames<-` renames the labels as well.

```
colnames<-,
rownames<-
```

## 7. Creating a hyperSpec Object, Data Import and Export

*hyperSpec* comes with filters for a variety of file formats. These are discussed in detail in a separate vignette accessible via `vignette ("file-io")`.

### 7.1. Creating a hyperSpec Object from Spectra Matrix and Wavelength Vector

If the data is in R's workspace, a *hyperSpec* object is created by:

```
> spc <- new ("hyperSpec", spc = spectra.matrix, wavelength = wavelength.vector, data = extra.data)
```

The most frequently needed arguments are:

```
spc          the spectra matrix
wavelength   the wavelength axis vector
data         the extra data (can already contain the spectra matrix in column $spc)
label        a list with the proper labels. Do not forget the wavelength axis label in $.wavelength
              and the spectral intensity axis label in $spc.
```

## 8. The Logbook

Slot `@log` of *hyperSpec* objects is intended to keep track of the history of the object. This logbook part of the output of the `summary`, and can also be retrieved by `logbook`.

```
> logbook (flu)
```

	short.description	long.description		date	user
1	scan.txt.PerkinElmer	rawdata/....	2011-01-15	18:24:51	cb@cb
2	\$<-	c, : [n....	2011-01-15	18:24:51	cb@cb
3	labels<-	c, c / (....	2011-01-15	18:24:51	cb@cb

New entries can be created manually by calling `logentry`:

```
> tmp <- logentry (flu, short = "test", long = "This could also be a list of parameters")
> logbook (tmp)
```

`logentry`,  
`logbook`

	short.description	long.description		date	user
1	scan.txt.PerkinElmer	rawdata/....	2011-01-15	18:24:51	cb@cb
2	\$<-	c, : [n....	2011-01-15	18:24:51	cb@cb
3	labels<-	c, c / (....	2011-01-15	18:24:51	cb@cb
4	test	This cou....	2011-02-13	22:10:51	cb@cb

In addition, *hyperSpec* by default logs automatically all changes to the object:

```
> tmp <- tmp [1:3]
> logbook (tmp)
```

	short.description	long.description		date	user
1	scan.txt.PerkinElmer	rawdata/....	2011-01-15	18:24:51	cb@cb
2	\$<-	c, : [n....	2011-01-15	18:24:51	cb@cb
3	labels<-	c, c / (....	2011-01-15	18:24:51	cb@cb
4	test	This cou....	2011-02-13	22:10:51	cb@cb
5	[]	i, j, .....	2011-02-13	22:10:51	cb@cb

The automatic logging mechanism can only log function calls and parameters (as opposed to the intention of the function call). *hyperSpec* functions that return a changed object allow to use more meaningful short descriptions: they are assigned via the argument *short*:

```
> tmp <- sweep (tmp, 2, mean, short = "centering")
> logbook (tmp)
```

	short.description	long.description		date	user
1	scan.txt.PerkinElmer	rawdata/....	2011-01-15	18:24:51	cb@cb
2	\$<-	c, : [n....	2011-01-15	18:24:51	cb@cb
3	labels<-	c, c / (....	2011-01-15	18:24:51	cb@cb
4	test	This cou....	2011-02-13	22:10:51	cb@cb
5	[]	i, j, .....	2011-02-13	22:10:51	cb@cb
6	centering	MARGIN, ....	2011-02-13	22:10:51	cb@cb

## 9. Combining and Decomposing hyperspec Objects

### 9.1. Binding Objects together

*hyperSpec* Objects can be bound together, either by columns (`cbind`) to append a new spectral range or by row (`rbind`) to append new spectra: `cbind` `rbind`

```
> dim (flu)
```

```
nrow ncol nwl
6 3 181
```

```
> dim (cbind (flu, flu))
```

```
nrow ncol nwl
6 3 362
```

```
> dim (rbind (flu, flu))
```

```
nrow ncol nwl
12 3 181
```



There is also a more general function, `bind`, taking the direction ("`r`" or "`c`") as first argument followed by the objects to bind either in separate arguments or in a list.

As usual for `rbind` and `cbind`, the objects that should be bound together must have the same rows and columns, respectively.

## 9.2. Binding Objects that do not Share the Same Extra Data and/or Wavelength Axis

`collapse` combines objects that should be bound together by row, but they do not share the columns and/or spectral range. The resulting object has all columns from all input objects, and all wavelengths from the input objects. If an input object does not have a particular column or wavelength, its value in the resulting object is `NA`.

`collapse`

The `barbituates` data is a list of 286 *hyperSpec* objects, each containing one mass spectrum. The spectra have between 4 and 101 data points each.

```
> barb <- collapse (barbituates)
> wl (barb) [1 : 25]

[1] 160.90 158.85 147.00 140.90 133.05 130.90 119.95 119.15 118.05 116.95 112.90 106.00 105.10
[14] 98.95 96.95 91.00 85.05 83.05 77.00 71.90 71.10 70.00 69.00 57.10 56.10
```

The resulting object does not have an ordered wavelength axis. This can be obtained in a second step:

```
> barb <- orderwl (barb)
> barb [[1:3, , min ~ min + 10i]]

      25.95 26.05 26.15 26.95 27.05 27.15 28.05 28.15 29.05 29.15 29.95
[1,]   NA   NA   NA   NA   562   NA   NA 11511  6146   NA   NA
[2,]   NA   NA   NA   NA   NA   NA   618 10151   NA  5040   NA   NA
[3,]   NA   NA   NA   NA   638   NA   NA 10722  5253   NA   NA
```

## 9.3. Binding Objects that do not Share the Same Spectra

`merge` adds a new spectral range (like `cbind`), but works also if spectra are missing in one of the objects. The arguments `by`, `by.x`, and `by.y` specify which columns should be used to decide which spectra are the same. The arguments `all`, `all.x`, and `all.y` determine whether spectra should be kept for the result set if they appear in only one of the objects. For details, see also the help on the base function `merge`.

`merge`

As an example, let's construct a version of the `chondro` data like being taken as two maps with different spectral ranges. In each data set, some spectra are missing.

```
> chondro.low <- sample (chondro [, , 600 ~ 1200], 700)
> nrow (chondro.low)

[1] 700

> chondro.high <- sample (chondro [, , 1400 ~ 1800], 700)
> nrow (chondro.high)

[1] 700
```

As all extra data columns are the same, no special declarations are needed for merging the data:

```
> chondro.merged <- merge (chondro.low, chondro.high)
> nrow (chondro.merged)

[1] 562
```

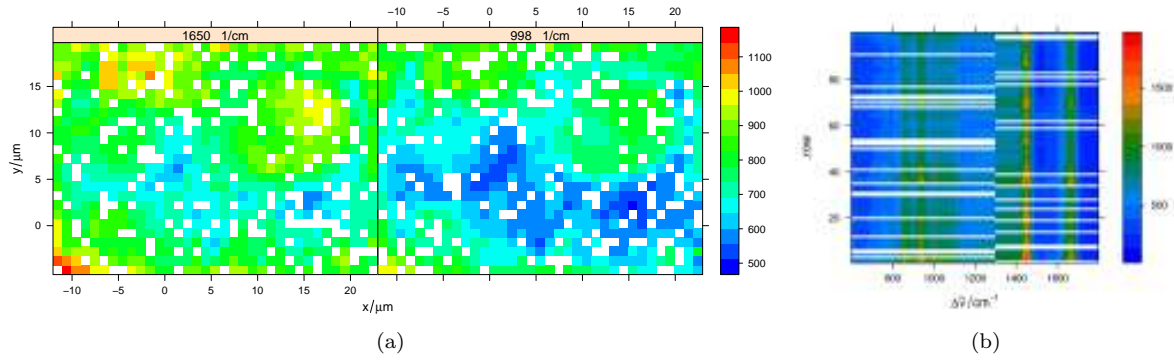


Figure 2: (a) For both spectral ranges some spectra are missing. (b) The missing parts of the spectra are filled with NA.

By default, the result consists of only those spectra, where *both* spectral ranges were available. To keep all spectra replacing missing parts by NA:

```
> chondro.merged <- merge (chondro.low, chondro.high, all = TRUE)
> nrow (chondro.merged)

[1] 838

> merged <- merge (chondro [1:7,, 610 ~ 620], chondro [5:10,, 615 ~ 625], all = TRUE)
> merged$.

```

	y	x	clusters	.nx	.ny	spc.610	spc.614	spc.618	spc.614	spc.618	spc.622	spc.626
1	-4.77	-11.55	matrix	1	NA	488.63	466.18	492.00	NA	NA	NA	NA
2	-4.77	-10.55	matrix	2	NA	489.48	465.05	490.53	NA	NA	NA	NA
3	-4.77	-9.55	matrix	3	NA	456.03	436.62	458.06	NA	NA	NA	NA
4	-4.77	-8.55	matrix	4	NA	464.82	444.85	470.02	NA	NA	NA	NA
5	-4.77	-7.55	matrix	5	1	428.66	410.80	433.12	410.80	433.12	461.19	397.38
6	-4.77	-6.55	matrix	6	2	426.07	407.86	431.21	407.86	431.21	458.15	394.18
7	-4.77	-5.55	lacuna	7	3	412.37	396.50	421.27	396.50	421.27	445.54	382.72
8	-4.77	-4.55	lacuna	NA	4	NA	NA	NA	381.95	406.25	429.67	368.46
9	-4.77	-3.55	lacuna	NA	5	NA	NA	NA	397.51	423.30	446.15	381.87
10	-4.77	-2.55	lacuna	NA	6	NA	NA	NA	377.39	402.23	424.19	362.43

If the spectra overlap, the result will have both data points. In the example here one could easily delete duplicate wavelengths. For real data, however, the duplicated wavelength will hardly ever contain the same values. The appropriate method to deal with this situation depends on the data at hand, but it will usually be some kind of spectral interpolation.

One possibility is removing duplicated wavelengths by using the mean intensity. This can conveniently be done by using `approx` using `method = "constant"`. For duplicated wavelengths, the intensities will be combined by the `tie` function. This already defaults to the mean, but we need `na.rm = TRUE`.

Thus, the function to calculate the new spectral intensities is

```
> approxfun <- function (y, wl, new.wl){
+   approx (wl, y, new.wl, method = "constant",
+         ties = function (x) mean (x, na.rm = TRUE)
+       )$y
+ }
```

which can be applied to the spectra:

```
> merged <- apply (merged, 1, approxfun,
+                 wl = wl (merged), new.wl = unique (wl (merged)),
+                 new.wavelength = "new.wl")
> merged$.
```

	y	x	clusters	.nx	.ny	spc.610	spc.614	spc.618	spc.622	spc.626
1	-4.77	-11.55	matrix	1	NA	488.6323....	466.1774....	492.0015....	NA	NA
2	-4.77	-10.55	matrix	2	NA	489.4758....	465.0506....	490.5328....	NA	NA
3	-4.77	-9.55	matrix	3	NA	456.0323....	436.6220....	458.0576....	NA	NA
4	-4.77	-8.55	matrix	4	NA	464.8207....	444.8485....	470.0171....	NA	NA
5	-4.77	-7.55	matrix	5	1	428.6619....	410.7955....	433.1227....	461.1903....	397.3773....
6	-4.77	-6.55	matrix	6	2	426.0734....	407.8569....	431.2144....	458.1502....	394.1775....
7	-4.77	-5.55	lacuna	7	3	412.3674....	396.5000....	421.2737....	445.5431....	382.7197....
8	-4.77	-4.55	lacuna	NA	4	NA	381.9504....	406.2470....	429.6728....	368.4599....
9	-4.77	-3.55	lacuna	NA	5	NA	397.5075....	423.3002....	446.1478....	381.8674....
10	-4.77	-2.55	lacuna	NA	6	NA	377.3917....	402.2348....	424.1901....	362.4296....

## 9.4. Matrix Multiplication

Two *hyperSpec* objects can be matrix multiplied by `%*%`. For an example, see the principal component analysis below (section 13.1 on page 29).

## 9.5. Decomposition

Matrix decompositions are common operations during chemometric data analysis. The results, e. g. of a principal component analysis are two matrices, the so-called scores and loadings. The results can have either the same number of rows as the spectra matrix they were calculated from (scores-like), or they have as many wavelengths as the spectra (loadings-like).

Both types of result objects can be “re-imported” into *hyperSpec* objects with function `decomposition`. A scores-like object retains all per-spectrum information (i. e. the extra data) while the spectra matrix and wavelength vector are replaced. A loadings-like object retains the wavelength information, while extra data is deleted (set to NA) unless the value is constant for all spectra.

A demonstration can be found in the principal component analysis example (section 13.1) on page 29.

## 10. Access to the data

The main functions to retrieve the data of a *hyperSpec* object are `[]` and `[] []`.


The difference between these functions is that `[]` returns a *hyperSpec* object, whereas the result of `[] []` is a *data.frame* if extra data columns were selected or otherwise the spectra matrix. Single extra data columns may be retrieved by `$`.

In order to change data, use `[]<-`, `[] []<-`, and `$<-` (see ).


### 10.1. Selecting and Deleting Spectra

The extraction function `[]` takes the spectra as first argument (For detailed help: see `?`[]``). It may be a vector giving the indices of the spectra to extract (select), a vector with negative indices indicating which spectra should be deleted, or a logical. Note that a matrix given to `[]` will be treated as a vector.


```
> plot (flu, col = "gray")
> plot (flu [1 : 3], add = TRUE)
```

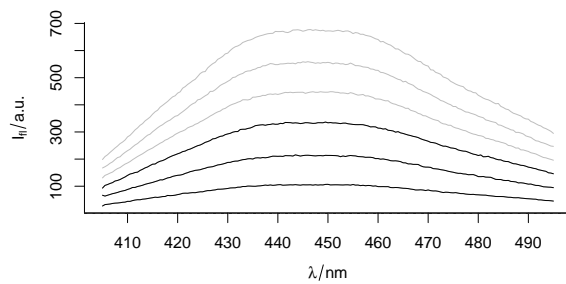

```
> plot (flu, col = "gray")
> plot (flu [-3], add = TRUE)
```

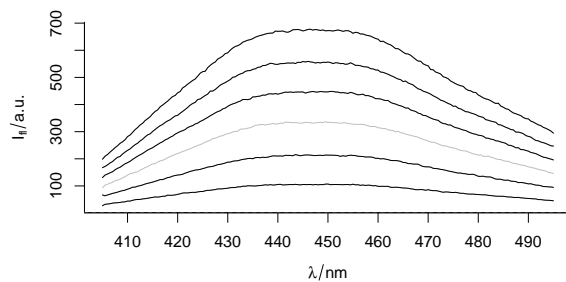

```
> plot (flu, col = "gray")
> plot (flu [flu$c > 0.2], add = TRUE)
```

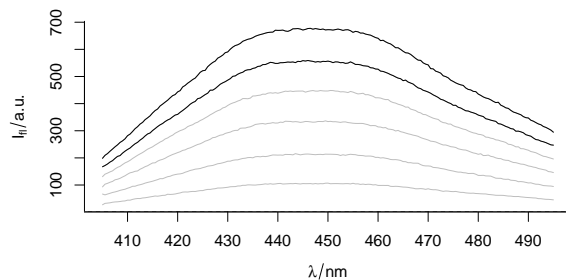

### 10.1.1. Random Samples

A random subset of spectra is conveniently selected by `sample` :

`sample`

```
> sample (chondro, 3)
```

```
hyperSpec object
 3 spectra
 4 data columns
300 data points / spectrum
wavelength: Delta * tilde(nu)/cm^-1 [numeric] 602 606 ... 1798
data: (3 rows x 4 columns)
 1. y: y/(mu * m) [numeric] -0.77 8.23 -4.77
 2. x: x/(mu * m) [numeric] 16.45 3.45 -11.55
 3. clusters: clusters [factor] matrix lacuna matrix
 4. spc: I / a.u. [matrix300] 292.42 293.72 ... 286.63
```

If appropriate indices into the spectra are needed instead, use `isample`:

`isample`

```
> isample (chondro, 3)
```

```
[1] 376 179 266
```

### 10.1.2. Sequences

Sequences of every  $n^{\text{th}}$  spectrum or the like can be retrieved with `seq`:

`seq`

```
> seq (chondro, length.out = 3, index = TRUE)
```

```
[1] 1 438 875
```

```
> seq (chondro, by = 100)
```

```
hyperSpec object
  9 spectra
  4 data columns
  300 data points / spectrum
wavelength: Delta * tilde(nu)/cm^-1 [numeric] 602 606 ... 1798
data: (9 rows x 4 columns)
  1. y: y/(mu * m) [numeric] -4.77 -2.77 ... 17.23
  2. x: x/(mu * m) [numeric] -11.55 18.45 ... 18.45
  3. clusters: clusters [factor] matrix matrix ... lacuna
  4. spc: I / a.u. [matrix300] 501.82 400.94 ... 124.64
```

Here, indices may be requested using `index = TRUE`.

### 10.2. Selecting Extra Data Columns

The second argument of the extraction functions `[]` and `[[[]]` specifies the (extra) data columns. They can be given like any column specification for a *data.frame*, i. e. numeric, logical, or by a vector of the column names:

They can be given like any column specification for a *data.frame*, i. e. numeric, logical, or by a vector of the column names:

```
> colnames (chondro)
```

```
[1] "y"          "x"          "clusters" "spc"
```

```
> chondro [[1 : 3, 1]]
```

```
      y
1 -4.77
2 -4.77
3 -4.77
```

```
> chondro [[1 : 3, -3]]
```

```
      y      x spc.602 spc.606 spc.610 spc.614 spc.618 spc.622 spc.626 spc.630 spc.634 spc.638
1 -4.77 -11.55 501.82 504.89 488.63 466.18 492.00 523.97 451.97 428.27 424.57 438.30
2 -4.77 -10.55 500.46 507.81 489.48 465.05 490.53 525.67 451.11 423.92 419.37 434.92
3 -4.77  -9.55 465.96 474.63 456.03 436.62 458.06 490.12 422.52 395.25 391.03 406.15
      spc.642 spc.646 spc.650 spc.654 spc.658 spc.662 spc.666 spc.670 spc.674 spc.678 spc.682 spc.686
1 481.30 460.99 447.44 450.24 450.94 453.38 458.20 466.69 469.25 473.57 482.20 488.52
2 480.81 462.10 443.69 446.64 447.99 451.27 456.25 464.99 466.85 471.44 480.55 487.78
3 448.59 431.83 414.82 416.81 416.96 419.79 423.57 432.40 435.08 441.93 448.94 454.50
      spc.690 spc.694 spc.698 spc.702 spc.706 spc.710 spc.714 spc.718 spc.722 spc.726 spc.730 spc.734
1 495.93 512.87 557.03 556.33 534.46 536.85 561.96 588.49 595.90 585.92 575.67 573.70
2 492.80 513.07 561.90 560.30 535.13 536.31 563.74 592.68 599.88 588.32 577.39 575.76
```

3	456.24	477.91	521.04	515.56	495.74	496.32	523.69	549.14	557.03	547.72	535.38	531.50
	spc.738	spc.742	spc.746	spc.750	spc.754	spc.758	spc.762	spc.766	spc.770	spc.774	spc.778	spc.782
1	583.16	607.77	625.35	652.34	689.63	718.75	711.28	690.89	670.42	658.17	647.64	639.97
2	587.49	614.27	630.52	658.00	695.37	728.06	716.36	696.24	677.02	663.74	650.89	642.73
3	540.20	565.21	583.46	603.63	637.19	666.98	657.48	642.32	626.72	611.11	600.15	593.50
	spc.786	spc.790	spc.794	spc.798	spc.802	spc.806	spc.810	spc.814	spc.818	spc.822	spc.826	spc.830
1	625.49	608.08	596.15	602.39	627.02	677.54	774.48	847.78	822.65	777.28	758.03	748.45
2	628.57	611.84	597.32	605.37	630.32	680.69	781.98	858.21	832.50	788.75	769.53	757.46
3	580.22	562.37	553.91	561.33	584.19	631.08	721.42	789.26	770.56	731.93	709.54	702.59
	spc.834	spc.838	spc.842	spc.846	spc.850	spc.854	spc.858	spc.862	spc.866	spc.870	spc.874	spc.878
1	752.22	790.32	871.64	1008.50	1177.04	1271.18	1235.81	1150.92	1109.33	1108.42	1141.15	1134.86
2	762.17	801.37	886.42	1026.68	1202.47	1298.98	1260.00	1172.42	1130.98	1128.85	1164.31	1158.62
3	706.89	741.03	819.82	947.87	1106.95	1201.04	1164.45	1086.59	1048.43	1044.20	1079.60	1073.95
	spc.882	spc.886	spc.890	spc.894	spc.898	spc.902	spc.906	spc.910	spc.914	spc.918	spc.922	spc.926
1	1070.75	1013.61	973.67	921.79	873.44	857.85	864.27	922.91	1040.13	1150.31	1171.55	1137.78
2	1089.56	1032.54	995.71	937.03	884.43	868.51	875.31	935.94	1055.43	1168.49	1192.43	1154.95
3	1010.01	964.35	934.66	878.11	824.49	810.28	814.78	869.52	978.28	1083.33	1105.26	1067.73
	spc.930	spc.934	spc.938	spc.942	spc.946	spc.950	spc.954	spc.958	spc.962	spc.966	spc.970	spc.974
1	1175.67	1312.52	1390.57	1342.27	1245.47	1180.60	1147.35	1133.25	1102.02	1077.46	1060.93	1021.80
2	1197.42	1341.38	1421.48	1371.01	1271.24	1205.06	1167.92	1151.16	1118.67	1095.37	1077.09	1040.31
3	1111.19	1239.52	1309.73	1268.45	1178.81	1115.80	1082.01	1070.06	1037.88	1019.11	997.35	964.71
	spc.978	spc.982	spc.986	spc.990	spc.994	spc.998	spc.1002	spc.1006	spc.1010	spc.1014	spc.1018	
1	966.40	900.98	872.01	877.06	874.93	982.55	1505.02	1243.86	959.96	911.68	872.39	
2	977.06	909.81	880.76	885.62	882.52	997.69	1553.27	1276.37	974.46	923.91	881.47	
3	909.89	850.79	825.91	823.99	821.05	930.51	1441.40	1177.69	904.44	857.50	821.46	
	spc.1022	spc.1026	spc.1030	spc.1034	spc.1038	spc.1042	spc.1046	spc.1050	spc.1054	spc.1058		
1	897.69	989.78	1129.23	1119.20	1025.25	989.58	979.85	973.84	983.36	1032.36		
2	906.69	1004.56	1149.62	1138.30	1037.90	1001.68	991.62	985.38	998.68	1057.71		
3	848.75	942.23	1075.50	1060.02	973.94	936.62	926.87	922.03	940.81	1006.63		
	spc.1062	spc.1066	spc.1070	spc.1074	spc.1078	spc.1082	spc.1086	spc.1090	spc.1094	spc.1098		
1	1123.79	1077.02	1001.11	981.62	989.00	997.05	1000.57	999.30	1008.65	999.79		
2	1161.68	1104.87	1016.74	992.43	1002.91	1010.31	1016.69	1017.98	1026.91	1013.99		
3	1120.71	1049.97	955.67	930.83	939.73	947.62	956.37	960.66	966.86	951.89		
	spc.1102	spc.1106	spc.1110	spc.1114	spc.1118	spc.1122	spc.1126	spc.1130	spc.1134	spc.1138		
1	981.23	935.15	873.40	822.73	807.88	848.08	930.19	899.08	761.12	702.87		
2	993.16	946.04	881.77	826.54	811.23	858.00	958.84	926.27	771.04	698.84		
3	929.49	887.57	827.85	772.72	759.32	810.09	916.48	889.30	729.66	657.33		
	spc.1142	spc.1146	spc.1150	spc.1154	spc.1158	spc.1162	spc.1166	spc.1170	spc.1174	spc.1178		
1	688.85	699.83	711.82	752.12	800.10	811.48	802.01	801.18	810.94	778.28		
2	682.10	694.61	705.31	749.96	799.95	810.09	797.99	799.33	812.78	775.15		
3	645.61	650.69	661.23	696.70	745.38	752.42	740.68	743.81	756.46	721.40		
	spc.1182	spc.1186	spc.1190	spc.1194	spc.1198	spc.1202	spc.1206	spc.1210	spc.1214	spc.1218		
1	734.79	713.11	726.36	757.89	793.61	850.44	899.87	872.93	833.37	855.55		
2	727.23	702.93	718.72	751.07	786.78	849.34	899.35	874.61	829.46	852.86		
3	677.04	652.16	667.27	702.84	733.76	785.85	834.72	813.09	772.95	789.94		
	spc.1222	spc.1226	spc.1230	spc.1234	spc.1238	spc.1242	spc.1246	spc.1250	spc.1254	spc.1258		
1	944.86	1075.75	1246.16	1398.67	1507.96	1563.22	1562.66	1513.41	1453.25	1413.05		
2	945.20	1083.59	1263.39	1419.18	1531.35	1587.42	1588.03	1537.70	1475.97	1433.96		
3	877.22	1001.95	1163.82	1308.93	1419.46	1462.29	1458.57	1415.37	1360.35	1324.38		
	spc.1262	spc.1266	spc.1270	spc.1274	spc.1278	spc.1282	spc.1286	spc.1290	spc.1294	spc.1298		
1	1404.31	1418.33	1405.56	1353.21	1291.22	1249.27	1208.26	1177.39	1230.13	1170.62		
2	1427.77	1445.69	1429.08	1371.05	1305.05	1264.35	1222.42	1195.63	1273.38	1209.68		
3	1320.97	1342.66	1328.15	1282.65	1216.06	1184.04	1147.58	1122.87	1224.06	1154.24		
	spc.1302	spc.1306	spc.1310	spc.1314	spc.1318	spc.1322	spc.1326	spc.1330	spc.1334	spc.1338		
1	1095.80	1087.78	1110.40	1142.35	1159.58	1123.49	1074.15	1053.12	1058.82	1080.48		
2	1113.95	1100.89	1122.56	1153.25	1171.15	1132.44	1083.39	1060.41	1068.13	1089.27		
3	1048.00	1025.97	1042.32	1067.15	1085.52	1047.97	1001.77	977.89	986.19	1007.82		
	spc.1342	spc.1346	spc.1350	spc.1354	spc.1358	spc.1362	spc.1366	spc.1370	spc.1374	spc.1378		
1	1079.82	1033.23	951.38	890.22	848.60	835.09	837.53	874.25	909.75	944.57		
2	1087.82	1035.05	952.08	887.65	841.71	826.79	830.62	870.13	903.96	939.22		
3	1008.66	955.74	881.38	822.10	777.17	764.62	770.80	808.37	844.40	873.10		
	spc.1382	spc.1386	spc.1390	spc.1394	spc.1398	spc.1402	spc.1406	spc.1410	spc.1414	spc.1418		
1	947.71	947.86	957.76	970.20	984.18	990.02	980.40	976.01	1007.97	1076.86		
2	938.48	938.33	946.12	961.70	974.15	979.99	970.29	968.31	999.53	1074.78		
3	869.52	867.73	874.67	884.85	895.82	902.04	893.01	890.76	918.58	995.42		

	spc.1422	spc.1426	spc.1430	spc.1434	spc.1438	spc.1442	spc.1446	spc.1450	spc.1454	spc.1458
1	1156.74	1199.57	1227.25	1367.10	1563.37	1642.74	1745.46	1801.44	1724.67	1630.80
2	1162.30	1210.52	1247.18	1408.17	1626.93	1698.05	1789.40	1850.11	1765.11	1672.56
3	1077.61	1117.95	1160.85	1320.58	1525.54	1572.53	1655.69	1707.09	1635.88	1564.66
	spc.1462	spc.1466	spc.1470	spc.1474	spc.1478	spc.1482	spc.1486	spc.1490	spc.1494	spc.1498
1	1535.44	1385.13	1184.59	978.49	796.08	665.70	563.59	503.36	461.91	440.43
2	1566.81	1407.78	1194.22	972.39	778.51	643.70	536.09	471.84	430.40	404.56
3	1462.17	1311.06	1113.78	902.97	722.46	595.63	498.79	435.84	401.15	376.79
	spc.1502	spc.1506	spc.1510	spc.1514	spc.1518	spc.1522	spc.1526	spc.1530	spc.1534	spc.1538
1	419.52	409.68	404.10	399.04	400.87	399.49	403.56	409.54	415.99	427.63
2	382.86	371.47	364.34	362.85	361.60	359.26	364.08	370.81	374.69	388.85
3	355.46	344.79	338.43	335.43	335.47	333.63	334.32	341.22	346.13	356.60
	spc.1542	spc.1546	spc.1550	spc.1554	spc.1558	spc.1562	spc.1566	spc.1570	spc.1574	spc.1578
1	443.99	470.50	493.26	497.57	484.15	473.95	467.82	466.12	466.70	470.97
2	405.99	433.79	457.91	460.94	443.97	433.85	429.54	426.30	428.09	431.84
3	373.42	398.94	419.96	421.96	407.62	397.15	394.34	393.09	393.67	397.23
	spc.1582	spc.1586	spc.1590	spc.1594	spc.1598	spc.1602	spc.1606	spc.1610	spc.1614	spc.1618
1	496.07	514.79	480.00	498.11	534.05	612.75	666.40	642.66	668.02	692.86
2	458.18	477.82	440.45	457.71	497.17	579.34	634.20	609.54	636.44	660.54
3	422.14	434.41	406.82	420.02	458.00	531.28	583.86	563.94	587.24	609.81
	spc.1622	spc.1626	spc.1630	spc.1634	spc.1638	spc.1642	spc.1646	spc.1650	spc.1654	spc.1658
1	728.91	799.08	881.58	966.82	1015.76	1056.97	1100.59	1154.02	1238.07	1332.49
2	694.97	767.17	853.22	938.43	990.41	1042.07	1085.39	1141.66	1227.04	1324.68
3	644.47	708.13	791.23	871.64	915.85	970.57	1014.74	1061.47	1138.61	1229.68
	spc.1662	spc.1666	spc.1670	spc.1674	spc.1678	spc.1682	spc.1686	spc.1690	spc.1694	spc.1698
1	1436.90	1513.97	1519.20	1396.05	1265.61	1164.23	1061.69	932.73	786.25	650.96
2	1429.39	1515.68	1526.92	1395.38	1253.24	1144.20	1034.00	899.96	748.38	607.60
3	1319.71	1400.43	1407.31	1296.09	1163.98	1057.79	959.04	830.98	693.73	564.39
	spc.1702	spc.1706	spc.1710	spc.1714	spc.1718	spc.1722	spc.1726	spc.1730	spc.1734	spc.1738
1	534.63	446.99	390.15	355.40	334.37	322.97	313.20	308.09	303.74	302.72
2	489.67	398.27	340.53	302.27	282.81	269.13	259.14	254.10	249.13	247.94
3	452.08	370.04	314.87	282.70	262.77	248.49	242.59	237.76	233.24	232.07
	spc.1742	spc.1746	spc.1750	spc.1754	spc.1758	spc.1762	spc.1766	spc.1770	spc.1774	spc.1778
1	297.49	296.82	295.90	291.05	290.78	287.34	290.81	286.37	286.36	284.07
2	241.96	243.22	241.29	235.20	234.71	232.18	234.36	230.91	230.17	226.88
3	227.71	225.45	225.13	219.18	219.37	216.38	218.99	216.25	215.78	212.53
	spc.1782	spc.1786	spc.1790	spc.1794	spc.1798					
1	287.07	285.02	285.01	286.12	286.63					
2	231.47	227.24	229.57	227.90	227.88					
3	215.52	213.42	214.86	212.83	212.80					

```
> chondro [[1 : 3, "x"]]
```

```

      x
1 -11.55
2 -10.55
3  -9.55

```

```
> chondro [[1 : 3, c(TRUE, FALSE, FALSE)]]
```

	y	spc.602	spc.606	spc.610	spc.614	spc.618	spc.622	spc.626	spc.630	spc.634	spc.638	spc.642
1	-4.77	501.82	504.89	488.63	466.18	492.00	523.97	451.97	428.27	424.57	438.30	481.30
2	-4.77	500.46	507.81	489.48	465.05	490.53	525.67	451.11	423.92	419.37	434.92	480.81
3	-4.77	465.96	474.63	456.03	436.62	458.06	490.12	422.52	395.25	391.03	406.15	448.59
	spc.646	spc.650	spc.654	spc.658	spc.662	spc.666	spc.670	spc.674	spc.678	spc.682	spc.686	spc.690
1	460.99	447.44	450.24	450.94	453.38	458.20	466.69	469.25	473.57	482.20	488.52	495.93
2	462.10	443.69	446.64	447.99	451.27	456.25	464.99	466.85	471.44	480.55	487.78	492.80
3	431.83	414.82	416.81	416.96	419.79	423.57	432.40	435.08	441.93	448.94	454.50	456.24
	spc.694	spc.698	spc.702	spc.706	spc.710	spc.714	spc.718	spc.722	spc.726	spc.730	spc.734	spc.738
1	512.87	557.03	556.33	534.46	536.85	561.96	588.49	595.90	585.92	575.67	573.70	583.16
2	513.07	561.90	560.30	535.13	536.31	563.74	592.68	599.88	588.32	577.39	575.76	587.49
3	477.91	521.04	515.56	495.74	496.32	523.69	549.14	557.03	547.72	535.38	531.50	540.20
	spc.742	spc.746	spc.750	spc.754	spc.758	spc.762	spc.766	spc.770	spc.774	spc.778	spc.782	spc.786
1	607.77	625.35	652.34	689.63	718.75	711.28	690.89	670.42	658.17	647.64	639.97	625.49
2	614.27	630.52	658.00	695.37	728.06	716.36	696.24	677.02	663.74	650.89	642.73	628.57
3	565.21	583.46	603.63	637.19	666.98	657.48	642.32	626.72	611.11	600.15	593.50	580.22

	spc.790	spc.794	spc.798	spc.802	spc.806	spc.810	spc.814	spc.818	spc.822	spc.826	spc.830	spc.834
1	608.08	596.15	602.39	627.02	677.54	774.48	847.78	822.65	777.28	758.03	748.45	752.22
2	611.84	597.32	605.37	630.32	680.69	781.98	858.21	832.50	788.75	769.53	757.46	762.17
3	562.37	553.91	561.33	584.19	631.08	721.42	789.26	770.56	731.93	709.54	702.59	706.89
	spc.838	spc.842	spc.846	spc.850	spc.854	spc.858	spc.862	spc.866	spc.870	spc.874	spc.878	spc.882
1	790.32	871.64	1008.50	1177.04	1271.18	1235.81	1150.92	1109.33	1108.42	1141.15	1134.86	1070.75
2	801.37	886.42	1026.68	1202.47	1298.98	1260.00	1172.42	1130.98	1128.85	1164.31	1158.62	1089.56
3	741.03	819.82	947.87	1106.95	1201.04	1164.45	1086.59	1048.43	1044.20	1079.60	1073.95	1010.01
	spc.886	spc.890	spc.894	spc.898	spc.902	spc.906	spc.910	spc.914	spc.918	spc.922	spc.926	spc.930
1	1013.61	973.67	921.79	873.44	857.85	864.27	922.91	1040.13	1150.31	1171.55	1137.78	1175.67
2	1032.54	995.71	937.03	884.43	868.51	875.31	935.94	1055.43	1168.49	1192.43	1154.95	1197.42
3	964.35	934.66	878.11	824.49	810.28	814.78	869.52	978.28	1083.33	1105.26	1067.73	1111.19
	spc.934	spc.938	spc.942	spc.946	spc.950	spc.954	spc.958	spc.962	spc.966	spc.970	spc.974	spc.978
1	1312.52	1390.57	1342.27	1245.47	1180.60	1147.35	1133.25	1102.02	1077.46	1060.93	1021.80	966.40
2	1341.38	1421.48	1371.01	1271.24	1205.06	1167.92	1151.16	1118.67	1095.37	1077.09	1040.31	977.06
3	1239.52	1309.73	1268.45	1178.81	1115.80	1082.01	1070.06	1037.88	1019.11	997.35	964.71	909.89
	spc.982	spc.986	spc.990	spc.994	spc.998	spc.1002	spc.1006	spc.1010	spc.1014	spc.1018	spc.1022	
1	900.98	872.01	877.06	874.93	982.55	1505.02	1243.86	959.96	911.68	872.39	897.69	
2	909.81	880.76	885.62	882.52	997.69	1553.27	1276.37	974.46	923.91	881.47	906.69	
3	850.79	825.91	823.99	821.05	930.51	1441.40	1177.69	904.44	857.50	821.46	848.75	
	spc.1026	spc.1030	spc.1034	spc.1038	spc.1042	spc.1046	spc.1050	spc.1054	spc.1058	spc.1062		
1	989.78	1129.23	1119.20	1025.25	989.58	979.85	973.84	983.36	1032.36	1123.79		
2	1004.56	1149.62	1138.30	1037.90	1001.68	991.62	985.38	998.68	1057.71	1161.68		
3	942.23	1075.50	1060.02	973.94	936.62	926.87	922.03	940.81	1006.63	1120.71		
	spc.1066	spc.1070	spc.1074	spc.1078	spc.1082	spc.1086	spc.1090	spc.1094	spc.1098	spc.1102		
1	1077.02	1001.11	981.62	989.00	997.05	1000.57	999.30	1008.65	999.79	981.23		
2	1104.87	1016.74	992.43	1002.91	1010.31	1016.69	1017.98	1026.91	1013.99	993.16		
3	1049.97	955.67	930.83	939.73	947.62	956.37	960.66	966.86	951.89	929.49		
	spc.1106	spc.1110	spc.1114	spc.1118	spc.1122	spc.1126	spc.1130	spc.1134	spc.1138	spc.1142		
1	935.15	873.40	822.73	807.88	848.08	930.19	899.08	761.12	702.87	688.85		
2	946.04	881.77	826.54	811.23	858.00	958.84	926.27	771.04	698.84	682.10		
3	887.57	827.85	772.72	759.32	810.09	916.48	889.30	729.66	657.33	645.61		
	spc.1146	spc.1150	spc.1154	spc.1158	spc.1162	spc.1166	spc.1170	spc.1174	spc.1178	spc.1182		
1	699.83	711.82	752.12	800.10	811.48	802.01	801.18	810.94	778.28	734.79		
2	694.61	705.31	749.96	799.95	810.09	797.99	799.33	812.78	775.15	727.23		
3	650.69	661.23	696.70	745.38	752.42	740.68	743.81	756.46	721.40	677.04		
	spc.1186	spc.1190	spc.1194	spc.1198	spc.1202	spc.1206	spc.1210	spc.1214	spc.1218	spc.1222		
1	713.11	726.36	757.89	793.61	850.44	899.87	872.93	833.37	855.55	944.86		
2	702.93	718.07	751.07	786.78	849.34	899.35	874.61	829.46	852.86	945.20		
3	652.16	667.27	702.84	733.76	785.85	834.72	813.09	772.95	789.94	877.22		
	spc.1226	spc.1230	spc.1234	spc.1238	spc.1242	spc.1246	spc.1250	spc.1254	spc.1258	spc.1262		
1	1075.75	1246.16	1398.67	1507.96	1563.22	1562.66	1513.41	1453.25	1413.05	1404.31		
2	1083.59	1263.39	1419.18	1531.35	1587.42	1588.03	1537.70	1475.97	1433.96	1427.77		
3	1001.95	1163.82	1308.93	1419.46	1462.29	1458.57	1415.37	1360.35	1324.38	1320.97		
	spc.1266	spc.1270	spc.1274	spc.1278	spc.1282	spc.1286	spc.1290	spc.1294	spc.1298	spc.1302		
1	1418.33	1405.56	1353.21	1291.22	1249.27	1208.26	1177.39	1230.13	1170.62	1095.80		
2	1445.69	1429.08	1371.05	1305.05	1264.35	1222.42	1195.63	1273.38	1209.68	1113.95		
3	1342.66	1328.15	1282.65	1216.06	1184.04	1147.58	1122.87	1224.06	1154.24	1048.00		
	spc.1306	spc.1310	spc.1314	spc.1318	spc.1322	spc.1326	spc.1330	spc.1334	spc.1338	spc.1342		
1	1087.78	1110.40	1142.35	1159.58	1123.49	1074.15	1053.12	1058.82	1080.48	1079.82		
2	1100.89	1122.56	1153.25	1171.15	1132.44	1083.39	1060.41	1068.13	1089.27	1087.82		
3	1025.97	1042.32	1067.15	1085.52	1047.97	1001.77	977.89	986.19	1007.82	1008.66		
	spc.1346	spc.1350	spc.1354	spc.1358	spc.1362	spc.1366	spc.1370	spc.1374	spc.1378	spc.1382		
1	1033.23	951.38	890.22	848.60	835.09	837.53	874.25	909.75	944.57	947.71		
2	1035.05	952.08	887.65	841.71	826.79	830.62	870.13	903.96	939.22	938.48		
3	955.74	881.38	822.10	777.17	764.62	770.80	808.37	844.40	873.10	869.52		
	spc.1386	spc.1390	spc.1394	spc.1398	spc.1402	spc.1406	spc.1410	spc.1414	spc.1418	spc.1422		
1	947.86	957.76	970.20	984.18	990.02	980.40	976.01	1007.97	1076.86	1156.74		
2	938.33	946.12	961.70	974.15	979.99	970.29	968.31	999.53	1074.78	1162.30		
3	867.73	874.67	884.85	895.82	902.04	893.01	890.76	918.58	995.42	1077.61		
	spc.1426	spc.1430	spc.1434	spc.1438	spc.1442	spc.1446	spc.1450	spc.1454	spc.1458	spc.1462		
1	1199.57	1227.25	1367.10	1563.37	1642.74	1745.46	1801.44	1724.67	1630.80	1535.44		
2	1210.52	1247.18	1408.17	1626.93	1698.05	1789.40	1850.11	1765.11	1672.56	1566.81		
3	1117.95	1160.85	1320.58	1525.54	1572.53	1655.69	1707.09	1635.88	1564.66	1462.17		
	spc.1466	spc.1470	spc.1474	spc.1478	spc.1482	spc.1486	spc.1490	spc.1494	spc.1498	spc.1502		



```

1 1385.13 1184.59 978.49 796.08 665.70 563.59 503.36 461.91 440.43 419.52
2 1407.78 1194.22 972.39 778.51 643.70 536.09 471.84 430.40 404.56 382.86
3 1311.06 1113.78 902.97 722.46 595.63 498.79 435.84 401.15 376.79 355.46
  spc.1506 spc.1510 spc.1514 spc.1518 spc.1522 spc.1526 spc.1530 spc.1534 spc.1538 spc.1542
1 409.68 404.10 399.04 400.87 399.49 403.56 409.54 415.99 427.63 443.99
2 371.47 364.34 362.85 361.60 359.26 364.08 370.81 374.69 388.85 405.99
3 344.79 338.43 335.43 335.47 333.63 334.32 341.22 346.13 356.60 373.42
  spc.1546 spc.1550 spc.1554 spc.1558 spc.1562 spc.1566 spc.1570 spc.1574 spc.1578 spc.1582
1 470.50 493.26 497.57 484.15 473.95 467.82 466.12 466.70 470.97 496.07
2 433.79 457.91 460.94 443.97 433.85 429.54 426.30 428.09 431.84 458.18
3 398.94 419.96 421.96 407.62 397.15 394.34 393.09 393.67 397.23 422.14
  spc.1586 spc.1590 spc.1594 spc.1598 spc.1602 spc.1606 spc.1610 spc.1614 spc.1618 spc.1622
1 514.79 480.00 498.11 534.05 612.75 666.40 642.66 668.02 692.86 728.91
2 477.82 440.45 457.71 497.17 579.34 634.20 609.54 636.44 660.54 694.97
3 434.41 406.82 420.02 458.00 531.28 583.86 563.94 587.24 609.81 644.47
  spc.1626 spc.1630 spc.1634 spc.1638 spc.1642 spc.1646 spc.1650 spc.1654 spc.1658 spc.1662
1 799.08 881.58 966.82 1015.76 1056.97 1100.59 1154.02 1238.07 1332.49 1436.90
2 767.17 853.22 938.43 990.41 1042.07 1085.39 1141.66 1227.04 1324.68 1429.39
3 708.13 791.23 871.64 915.85 970.57 1014.74 1061.47 1138.61 1229.68 1319.71
  spc.1666 spc.1670 spc.1674 spc.1678 spc.1682 spc.1686 spc.1690 spc.1694 spc.1698 spc.1702
1 1513.97 1519.20 1396.05 1265.61 1164.23 1061.69 932.73 786.25 650.96 534.63
2 1515.68 1526.92 1395.38 1253.24 1144.20 1034.00 899.96 748.38 607.60 489.67
3 1400.43 1407.31 1296.09 1163.98 1057.79 959.04 830.98 693.73 564.39 452.08
  spc.1706 spc.1710 spc.1714 spc.1718 spc.1722 spc.1726 spc.1730 spc.1734 spc.1738 spc.1742
1 446.99 390.15 355.40 334.37 322.97 313.20 308.09 303.74 302.72 297.49
2 398.27 340.53 302.27 282.81 269.13 259.14 254.10 249.13 247.94 241.96
3 370.04 314.87 282.70 262.77 248.49 242.59 237.76 233.24 232.07 227.71
  spc.1746 spc.1750 spc.1754 spc.1758 spc.1762 spc.1766 spc.1770 spc.1774 spc.1778 spc.1782
1 296.82 295.90 291.05 290.78 287.34 290.81 286.37 286.36 284.07 287.07
2 243.22 241.29 235.20 234.71 232.18 234.36 230.91 230.17 226.88 231.47
3 225.45 225.13 219.18 219.37 216.38 218.99 216.25 215.78 212.53 215.52
  spc.1786 spc.1790 spc.1794 spc.1798
1 285.02 285.01 286.12 286.63
2 227.24 229.57 227.90 227.88
3 213.42 214.86 212.83 212.80

```

To select one column, the `$` operator is more convenient:

`$`

```
> flu$c
```

```
[1] 0.05 0.10 0.15 0.20 0.25 0.30
```

*hyperSpec* supports command line completion for the `$` operator.

The extra data may also be set this way:

`$<-`


```
> flu$n <- list (1 : 6, label = "sample no.")
```

This function will append new columns, if necessary.

### 10.3. Selecting Wavelength Ranges


Wavelength ranges can easily be selected using `[]`'s third argument:

```
> plot (paracetamol [, , 2800 ~ 3200])
```

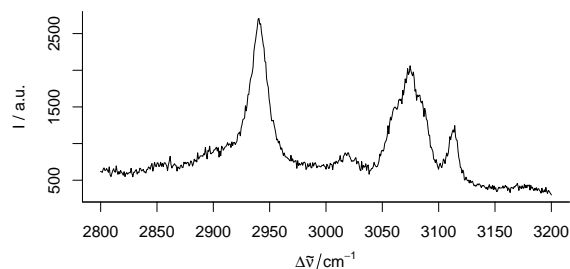

By default, the values given are treated as wavelengths, if they are indices into the columns of the spectra matrix, use `wl.index = TRUE`:

```
> plot (paracetamol [, 2800 : 3200, wl.index = TRUE])
```

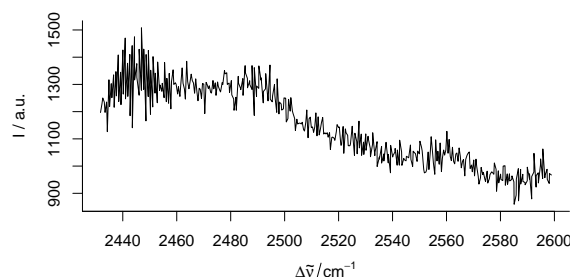


Section 10.4.1 (p. 19) details into the different possibilities of specifying wavelengths.

#### 10.4. Deleting Wavelength Ranges

Deleting wavelength ranges may be accomplished using negative index vectors together with `wl.index = TRUE`.


```
> plot (paracetamol [, -(500 : 1000), wl.index = TRUE])
```

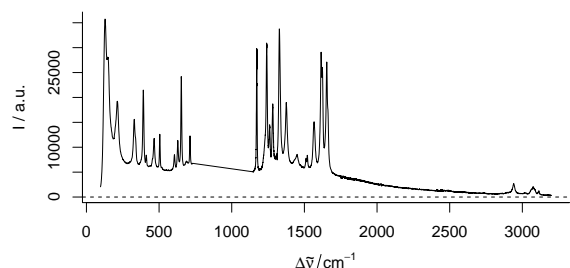

However, this mechanism works only if the proper indices are known.

If the range to be cut out is rather known in the units of the wavelength axis, it is easier to select the remainder of the spectrum instead. To delete the spectral range from 1750 to 2800  $\text{cm}^{-1}$  of the paracetamol spectrum one can thus use:

```
> plot (paracetamol [, c (min ~ 1750, 2800 ~ max)])
```

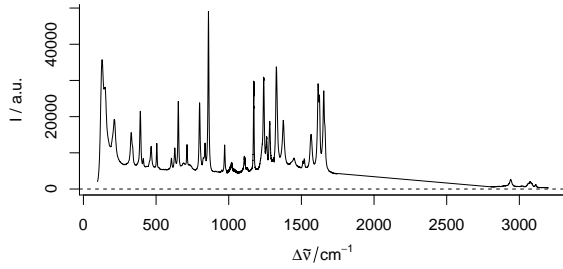

(It is possible to produce a plot of this data where the cut range is not bridged by a line and the wavelength axis is cut in order to save space. For details see the “plotting” vignette).

#### 10.4.1. Converting Wavelengths to Indices and vice versa

Spectra in *hyperSpec* have always discretized wavelength axes, they are stored in a matrix with column corresponding to one wavelength. *hyperSpec* provides two conversion functions: `i2w1` returns the wavelength corresponding to the given indices and `w12i` calculates index vectors from wavelengths.

`w12i i2w1`

If the wavelengths are given as a numeric vector, they are each converted to the corresponding wavelength. In addition there is a more sophisticated possibility of specifying wavelength ranges using a *formula*. The basic syntax is `start ~ end`. This yields a vector *index of start : index of end*.

The result of the formula conversion differs from the numeric vector conversion in three ways:

- The colon operator for constructing vectors accepts only integer numbers, the tilde (for formulas) does not have this restriction.
- If the vector does not take into account the spectral resolution, one may get only every  $n^{\text{th}}$  point or repetitions of the same index:

```
> w12i (flu, 405 : 410)
[1] 1 3 5 7 9 11

> w12i (flu, 405 ~ 410)
[1] 1 2 3 4 5 6 7 8 9 10 11

> w12i (chondro, 1000 : 1010)
[1] 100 101 101 101 101 101 102 102 102 102 103 103

> w12i (chondro, 1000 ~ 1010)
[1] 100 101 102 103
```

- If the object’s wavelength axis is not ordered, the formula approach will give weird results. In that (probably rare) case, use `orderw1` first to obtain an object with ordered wavelength axis.

*start* and *end* may contain the special variables `min` and `max` that correspond to the lowest and highest wavelengths of the object:

```
> w12i (flu, min ~ 410)
[1] 1 2 3 4 5 6 7 8 9 10 11
```

Often, specifications like *wavelength ± n data points* are needed. They can be given using complex numbers in the formula. The imaginary part is added to the index calculated from the wavelength in the real part:

```
> wl2i (flu, 450 - 2i ~ 450 + 2i)
```

```
[1] 89 90 91 92 93
```

```
> wl2i (flu, max - 2i ~ max)
```

```
[1] 179 180 181
```

To specify several wavelength ranges, use a list containing the formulas and vectors<sup>1</sup>:

```
> wl2i (flu, 450 - 2i ~ 450 + 2i)
```

```
[1] 89 90 91 92 93
```

```
> wl2i (flu, c (min ~ 406.5, max - 2i ~ max))
```

```
[1] 1 2 3 4 179 180 181
```

This mechanism also works for the wavelength arguments of `[]`, `[[ ]]`, and `plotspc`.

#### 10.4.2. Changing the Wavelength Axis

Sometimes wavelength axes need to be transformed, e. g. converting from wavelengths to frequencies. In this case, retrieve the wavelength axis vector with `wl`, convert each value of the resulting vector and assign the result with `wl<-`. Also the label of the wavelength axis may need to be adjusted.

`wl, wl<-`

As an example, convert the wavelength axis of `laser` to frequencies. As the wavelengths are in nanometers, and the frequencies are easiest expressed in terahertz, an additional conversion factor of 1000 is needed:

```
> laser
```

```
hyperSpec object
  84 spectra
  2 data columns
  36 data points / spectrum
wavelength: lambda/nm [numeric] 404.58 404.62 ... 405.82
data: (84 rows x 2 columns)
  1. t: t / s [numeric] 0 2 ... 5722
  2. spc: I / a.u. [matrix36] 164.65 179.72 ... 112.09
```

```
> wavelengths <- wl (laser)
```

```
> frequencies <- 2.998e8 / wavelengths / 1000
```

```
> wl (laser) <- frequencies
```

```
> labels (laser, ".wavelength") <- "f / THz"
```

```
> laser
```

```
hyperSpec object
  84 spectra
  2 data columns
  36 data points / spectrum
wavelength: f / THz [numeric] 741.01 740.95 ... 738.76
data: (84 rows x 2 columns)
  1. t: t / s [numeric] 0 2 ... 5722
  2. spc: I / a.u. [matrix36] 164.65 179.72 ... 112.09
```

```
> rm (laser)
```

---

<sup>1</sup>Formulas are combined to a list by `c`.

There are other possibilities of invoking `wl<-` including the new label, e.g.

```
> wl (laser, "f / THz") <- frequencies
and
> wl (laser) <- list (wl = frequencies, label = "f / THz")
see ?`wl<-` for more information.
```

### 10.4.3. Ordering the Wavelength Axis

If the wavelength axis of an object needs reordering (e.g. after `collapse`), `orderwl` can be used:

`orderwl`

```
> barb <- collapse (barbituates [1 : 3])
> wl (barb)

[1] 160.90 158.85 147.00 140.90 133.05 130.90 119.95 119.15 118.05 116.95 112.90 106.00 105.10
[14] 98.95 96.95 91.00 85.05 83.05 77.00 71.90 71.10 70.00 69.00 57.10 56.10 55.00
[27] 43.85 43.05 41.10 40.10 39.00 32.15 31.15 30.05 29.05 28.15 27.05 132.95 131.00
[40] 120.05 119.05 117.95 113.00 105.90 82.95 72.00 69.10 56.00 44.05 40.00 30.15 28.05
[53] 27.15 84.15 68.90 55.10 43.95

> barb <- orderwl (barb)
> wl (barb)

[1] 27.05 27.15 28.05 28.15 29.05 30.05 30.15 31.15 32.15 39.00 40.00 40.10 41.10
[14] 43.05 43.85 43.95 44.05 55.00 55.10 56.00 56.10 57.10 68.90 69.00 69.10 70.00
[27] 71.10 71.90 72.00 77.00 82.95 83.05 84.15 85.05 91.00 96.95 98.95 105.10 105.90
[40] 106.00 112.90 113.00 116.95 117.95 118.05 119.05 119.15 119.95 120.05 130.90 131.00 132.95
[53] 133.05 140.90 147.00 158.85 160.90
```

### 10.5. More on the Square-Bracket Operators for Replacing Values

`[[ ]]` also accepts index matrices of size  $n \times 2$ . In this case, a vector of values from the spectra matrix is returned.

```
> indexmatrix <- matrix (c (1 : 3, 1 : 3), ncol = 2)
> indexmatrix

      [,1] [,2]
[1,]    1    1
[2,]    2    2
[3,]    3    3

> chondro [[indexmatrix, wl.index = TRUE]]
[1] 501.82 507.81 456.03

> diag (chondro [[1 : 3, , min ~ min + 2i]])
[1] 501.82 507.81 456.03

[[ ]]<- also accepts index matrices of size  $n \times 2$ .

> indexmatrix <- matrix (c (1 : 3, 1 : 3), ncol = 2)
> indexmatrix

      [,1] [,2]
[1,]    1    1
[2,]    2    2
[3,]    3    3

> chondro [[indexmatrix, wl.index = TRUE]]
[1] 501.82 507.81 456.03

> diag (chondro [[1 : 3, , min ~ min + 2i]])
[1] 501.82 507.81 456.03
```

## 10.6. Fast Access to Parts of the hyperSpec Object

[[]] \$. \$..

*hyperSpec* comes with three abbreviation functions for easy access to the data:

```
x [[]] returns the spectra matrix (x$spc).
x [[i, , 1]] the cut spectra matrix is returned if wavelengths are specified in l.
x [[i, j, 1]] If data columns are selected (second index), the result is a data.frame.
x [[i, , 1]] <- Also, parts of the spectra matrix can be set (only indices for spectra and wavelength are allowed for this function).
x [i, j] <- sets parts of x@data.
x $. returns the complete data.frame x@data, with the spectra in column $spc.
x $.. returns the extra data (x@data without x$spc).
x $.. <- sets the extra data (x@data without x$spc). However, the columns must match exactly in this case.
```

## 10.7. Conversion to Long-Format data.frame

Some functions need the data being an *unstacked* or *long-format data.frame*. `as.long.df` is the appropriate conversion function. as.long.df

## 11. Plotting

*hyperSpec* offers a variety of possibilities to plot spectra, spectral maps, the spectra matrix, time series, depth profiles, etc.. This all is discussed in a separate document: see `vignette ("plotting")`.

## 12. Spectral (Pre)processing

### 12.1. Cutting the Spectral Range

[] [[]]

The extraction functions `[]` and `[[]]` can be used to cut the spectra: Their third argument takes wavelength specifications as discussed above and also logicals (i.e. vectors specifying with TRUE/FALSE for each column of `$spc` whether it should be included or not).

`[]` returns a *hyperSpec* object, `[[]]` the spectra matrix `$spc` (or the *data.frame* `@data` if in addition data columns were specified) only.

```
> flu [, , min ~ 408.5]

hyperSpec object
  6 spectra
  4 data columns
  8 data points / spectrum
wavelength: lambda/nm [numeric] 405.0 405.5 ... 408.5
data: (6 rows x 4 columns)
  1. file: [factor] rawdata/flu1.txt rawdata/flu2.txt ... rawdata/flu6.txt
  2. spc: I[fl]/a.u. [AsIs matrix x 8] 27.150 66.801 ... 256.89
  3. c: c / (mg / l) [numeric] 0.05 0.10 ... 0.3
  4. n: sample no. [integer] 1 2 ... 6

> flu [[, , c (min ~ min + 2i, max - 2i ~ max)]]
```

	405	405.5	406	494	494.5	495
[1,]	27.150	32.345	33.379	47.163	46.412	45.256
[2,]	66.801	63.715	66.712	96.602	96.206	94.610
[3,]	93.144	103.068	106.194	149.539	148.527	145.793
[4,]	130.664	139.998	143.798	201.484	198.867	195.867
[5,]	167.267	171.898	177.471	252.066	248.067	246.952
[6,]	198.430	209.458	215.785	307.519	302.325	294.649

## 12.2. Shifting Spectra

Sometimes, spectra need to be aligned along the spectral axis.

In general, two options are available for shifting spectra along the wavelength axis.

1. The wavelength axis can be shifted, while the intensities stay unaffected.
2. the spectra are interpolated onto a new wavelength axis, while the nominal wavelengths stay.

The first method is very straightforward:

```
> tmp <- chondro
> wl (tmp) <- wl (tmp) - 10
```

but it cannot be used if each spectrum (or groups of spectra) are shifted individually.

In that case, interpolation is needed. R offers many possibilities to interpolate (e.g. `approx` for constant / linear approximation, `spline` for spline interpolation, `loess` can be used to obtain smoothed approximations, etc.). The appropriate interpolation strategy will depend on the spectra, and *hyperSpec* therefore leaves it up to the user to select a sensible interpolation function.

As an example, we will use natural splines to do the interpolation. It is convenient to set it up as a function:

```
> interpolate <- function (spc, shift, wl){
+   spline (wl + shift, spc, xout = wl, method = "natural")$y
+ }
```

This function can now be applied to a set of spectra:

```
> tmp <- apply (chondro, 1, interpolate, shift = -10, wl = wl (chondro))
```

If different spectra need to be offset by different shift, use a loop<sup>2</sup>

```
> shifts <- rnorm (nrow (chondro))
> tmp <- chondro [[]
> for (i in seq_len (nrow (chondro)))
+   tmp [i, ] <- interpolate (tmp [i, ], shifts [i], wl = wl (chondro))
> chondro [[] <- tmp
```

### 12.2.1. Calculating the Shift

Often, the shift in the spectra is determined by aligning a particular signal. This strategy works best with spectrally oversampled data that allows accurate determination of the signal position.

For the `chondro` data, let's use the maximum of the phenylalanine band between 990 and 1020  $\text{cm}^{-1}$ . As just the very maximum is too coarse, we'll use the maximum of a square polynomial fitted to the maximum and its two neighbours.

<sup>2</sup>`sweep` cannot be used here, and while there is the possibility to use `sapply` or `mapply`, they are not faster than the for loop in this case. Make sure to work on a copy of the spectra matrix, as that is much faster than row-wise extracting and changing the spectra by `[]` and `[]<-`.

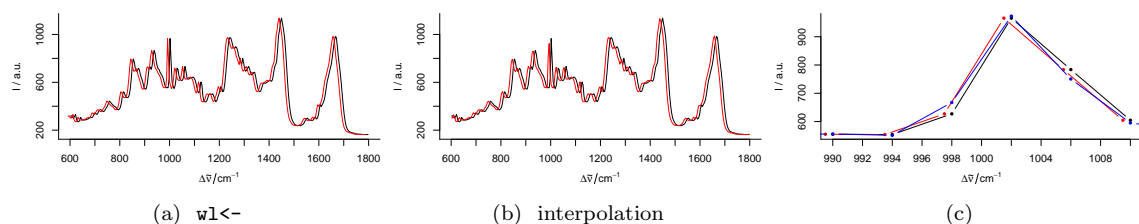


Figure 3: Shifting the Spectra along the Wavelength Axis. (a) Changing the wavelength values. (b) Interpolation. (c) Detail view of the phenylalanine band: shifting by `wl<-` (red) does not affect the intensities, while the spectrum is slightly changed by interpolations (blue).

```
> find.max <- function (y, x){
+   pos <- which.max (y) + (-1:1)
+   X <- x [pos] - x [pos [2]]
+   Y <- y [pos] - y [pos [2]]
+
+   X <- cbind (1, X, X^2)
+   coef <- qr.solve (X, Y)
+
+   - coef [2] / coef [3] / 2 + x [pos [2]]
+ }
> bandpos <- apply (chondro [[, 990 ~ 1020]], 1, find.max, wl (chondro [, 990 ~ 1020]))
> refpos <- find.max (colMeans (chondro [[, 990 ~ 1020]]), wl (chondro [, 990 ~ 1020]))
> shift1 <- refpos - bandpos
```

A second possibility is to optimize the shift. For this strategy, the spectra must be sufficiently similar, while low spectral resolution is compensated by using larger spectral windows.

```
> chondro <- chondro - spc.fit.poly.below (chondro [, min+3i ~ max - 3i], chondro)

Fitting with npts.min = 15

> chondro <- sweep (chondro, 1, rowMeans (chondro [[]], na.rm = TRUE), "/")

> targetfn <- function (shift, wl, spc, targetspc){
+   error <- spline (wl + shift, spc, xout = wl)$y - targetspc
+   sum (error^2)
+ }
> shift2 <- numeric (nrow (chondro))
> tmp <- chondro [[]]
> target <- colMeans (chondro [[]])
> for (i in 1 : nrow (chondro))
+   shift2 [i] <- unlist (optimize (targetfn, interval = c (-5, 5), wl = chondro@wavelength,
+                                   spc = tmp[i,], targetspc = target)$minimum)
```

Figure 4 shows that the second correction method works better for the chondrocyte data. This was expected, as the spectra are hardly or not oversampled, but are very similar to each other.

### 12.3. Smoothing Interpolation

Spectra acquired by grating instruments are frequently interpolated onto a new wavelength axis, e.g. because the unequal data point spacing should be removed. Also, the spectra can be smoothed: reducing the spectral resolution allows to increase the signal to noise ratio. For chemometric data

spc.bin  
spc.loess



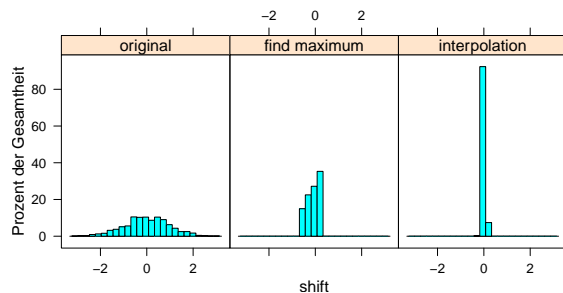


Figure 4: The shifts used to disturb the chondrocyte data (original), and the remaining shift after correction with the two methods discussed here.

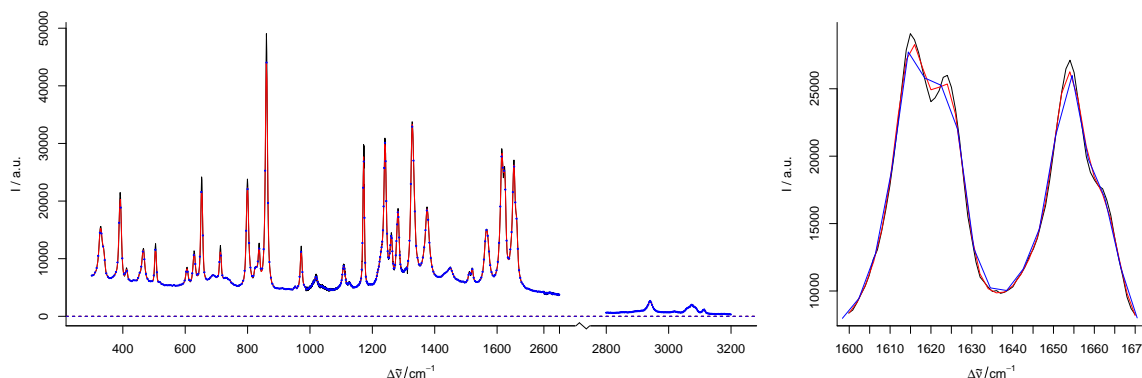


Figure 5: Smoothing interpolation by `spc.loess` with new data point spacing of  $2 \text{ cm}^{-1}$  (red) and `spc.bin` (blue). The magnification on the right shows how interpolation may cause a loss in signal height.

analysis reducing the number of data points per spectrum may be crucial as it reduces the dimensionality of the data.

*hyperSpec* provides two functions to do so: `spc.bin` and `spc.loess`.

`spc.bin` bins the spectral axis by averaging every *by* data points.

```
> plot (paracetamol, wl.range = c (300 ~ 1800, 2800 ~ max), xoffset = 850)
> p <- spc.loess (paracetamol, c(seq (300, 1800, 2), seq (2850, 3150, 2)))
> plot (p, wl.range = c (300 ~ 1800, 2800 ~ max), xoffset = 850, col = "red", add = TRUE)
> b <- spc.bin (paracetamol, 4)
> plot (b, wl.range = c (300 ~ 1800, 2800 ~ max), xoffset = 850,
+       lines.args = list (pch = 20, cex = .3, type = "p"), col = "blue", add = TRUE)
```

`spc.loess` applies R's `loess` function for spectral interpolation. Figure 5 shows the result of interpolating from 300 to 1800 and 2850 to  $3150 \text{ cm}^{-1}$  with  $2 \text{ cm}^{-1}$  data point distance. This corresponds to a spectral resolution of about  $4 \text{ cm}^{-1}$ , and the decrease in spectral resolution can be seen at the sharp bands where the maxima are not reached (due to the fact that the interpolation wavelength axis does not necessarily hit the maxima). The original spectrum had 4064 data points with unequal data point spacing (between 0 and  $1.4 \text{ cm}^{-1}$ ). The interpolated spectrum has 902 data points.

## 12.4. Background Correction

sweep

To subtract a background spectrum of each of the spectra in an object, use `sweep (spectra, 2, background.spectrum, "-")`.

## 12.5. Offset Correction

apply sweep

Calculate the offsets and sweep them off the spectra:

```
> offsets <- apply (chondro, 1, min)
> chondro.offset.corrected <- sweep (chondro, 1, offsets, "-")
```

If the offset is calculated by a function, as here with the `min`, *hyperSpec*'s `sweep` method offers a shortcut: `sweep`'s *STATS* argument may be the function instead of a numeric vector:

```
> chondro.offset.corrected <- sweep (chondro, 1, min, "-")
```

## 12.6. Baseline Correction

*hyperSpec* comes with two functions to fit polynomial baselines.

spc.fit.poly  
spc.fit.poly.below

`spc.fit.poly` fits a polynomial baseline of the given order. A least-squares fit is done so that the function may be used on rather noisy spectra. However, the user must supply an object that is cut appropriately. Particularly, the supplied wavelength ranges are not weighted.

`spc.fit.poly.below` tries to find appropriate support points for the baseline iteratively.

Both functions return a *hyperSpec* object containing the fitted baselines. They need to be subtracted afterwards:

```
> bl <- spc.fit.poly.below (chondro)
Fitting with npts.min = 15
> chondro <- chondro - bl
```

For details, see vignette (baselinebelow).

## 12.7. Intensity Calibration

### 12.7.1. Correcting by a constant, e. g. Readout Bias

CCD cameras often operate with a bias, causing a constant value for each pixel. Such a constant can be immediately subtracted:

```
spectra - constant
```

### 12.7.2. Correcting Wavelength Dependence

sweep

For each of the wavelengths the same correction needs to be applied to all spectra.

1. There might be wavelength dependent offsets (background or dark spectra). They are subtracted:  
`sweep (spectra, 2, offset.spectrum, "-")`
2. A multiplicative dependency such as a CCD's photon efficiency:  
`sweep (spectra, 2, photon.efficiency, "/")`

### 12.7.3. Spectra Dependent Correction

sweep

If the correction depends on the spectra (e.g. due to inhomogeneous illumination while collecting imaging data, differing optical path length, etc.), the *MARGIN* of the **sweep** function needs to be 1 or SPC:

1. Pixel dependent offsets are subtracted:  
`sweep (spectra, SPC, pixel.offsets, "-")`
2. A multiplicative dependency:  
`sweep (spectra, SPC, illumination.factors, "*")`

### 12.8. Normalization

apply sweep

Again, **sweep** is the function of choice. E.g. for area normalization, use:

```
> chondro <- sweep (chondro, 1, mean, "/")
```

(using the mean instead of the sum results in conveniently scaled spectra with intensities around 1.)

If the calculation of the normalization factors is more elaborate, use a two step procedure:

1. Calculate appropriate normalization factors  
You may calculate the factors using only a certain wavelength range, thereby normalizing on a particular band or peak.
2. Again, sweep the factor off the spectra:  
`normalized <- sweep (spectra, 1, factors, "*")`

```
> factors <- 1 / apply (chondro [, , 1600 ~ 1700], 1, mean)
```

```
> chondro <- sweep (chondro, 1, factors, "*")
```

For minimum-maximum-normalization, first do an offset- or baseline correction, then normalize using **max**.

### 12.9. Centering the Data


apply sweep

Centering means that the mean spectrum is subtracted from each of the spectra. Many data analysis techniques, like principal component analysis, partial least squares, etc., work much better on centered data.

However, from a spectroscopic point of view it depends on the particular data set whether centering does make sense or not.

To centre the **flu** data set, use:


```
> flu.centered <- sweep (flu, 2, mean, "-")  
> plot (flu.centered)
```

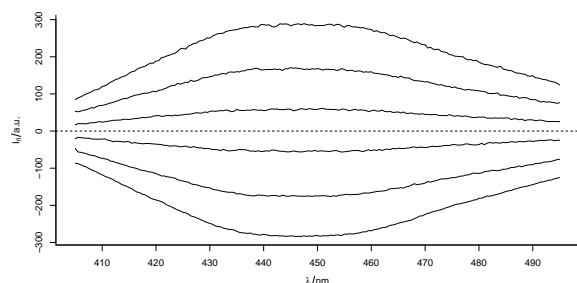

On the other hand, the `chondro` data set consists of Raman spectra, so the spectroscopic interpretation of centering is getting rid of the the average chemical composition of the sample. But: what is the meaning of the “average spectrum” of an inhomogeneous sample? In this case it may be better to subtract the minimum spectrum (which will hopefully have almost the same benefit on the data analysis) as it is the spectrum of that chemical composition that is underlying the whole sample.

One more point to consider is that the actual minimum spectrum will pick up (negative) noise. In order to avoid that, using e.g. the 5<sup>th</sup> percentile spectrum is more suitable:

```
> perc.5th <- apply (chondro, 2, quantile, 0.05)
> chondro <- sweep (chondro, 2, perc.5th, "-")
> plot (chondro, "spcprct15")
```

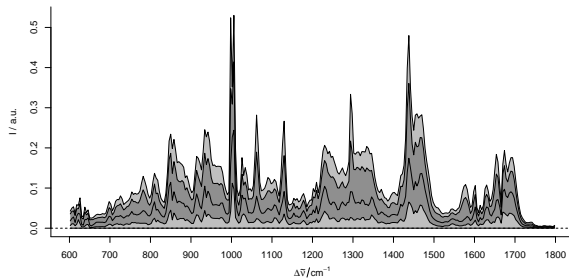

## 12.10. Variance Scaling

Variance scaling is often used in multivariate analysis to adjust the influence and scaling of the variates (that are typically different physical values). However, spectra already do have the same scale of the same physical value. Thus one has to trade off the the expected numeric benefit with the fact that wavelengths with low signal will contain exploded noise after variance scaling.

Again, `sweep` may be used:

```
> scaled.chondro <- sweep (chondro, 2, var, "/")
```

Alternatively, R provides a function `scale` which works on matrices:

```
> scaled.chondro <- chondro
> scaled.chondro [[]] <- scale (scaled.chondro [[]])
```

apply sweep  
scale

## 12.11. Multiplicative Scatter Correction (MSC)

MSC can be done using `msc` from package *pls*[1]. It operates on the spectra matrix:

```
> library (pls)
> chondro.msc <- chondro
> chondro.msc [[]] <- msc (chondro [[]])
```

pls::msc

## 12.12. Spectral Arithmetic

Basic mathematical functions are defined for *hyperSpec* objects. You may convert spectra:  
`absorbance.spectra = - log10 (transmission.spectra)`

In this case, do not forget to adapt the label:

```
> labels (absorbance.spectra)$spc <- "A"
```

+ - \* / ^ log  
log10

labels

Be careful: R's `log` function calculates the natural logarithm if no base is given.

The basic arithmetic operators work element-wise in R. Thus they all need either a scalar, or a matrix (or *hyperSpec* object) of the correct size.

Matrix multiplication is done by `%*%`, again each of the operands may be a matrix or a *hyperSpec* object, and must have the correct dimensions.

## 13. Data Analysis

### 13.1. Data Analysis Methods using a `data.frame`

#### e. g. Principal Component Analysis with `prcomp`

`$.`

The `$.` notation is handy, if a data analysis function expects a *data.frame*. The column names can then be used in the formula:

```
> pca <- prcomp (~ spc, data = chondro$, center = FALSE)
```

Results of such a decomposition can be put again into *hyperSpec* objects. This allows to plot e.g. the loading like spectra, or score maps, see figure 6.

```
> scores <- decomposition (chondro, pca$x, label.wavelength = "PC",
+                           label.spc = "score / a.u.")
> scores
```

```
hyperSpec object
  875 spectra
   4 data columns
  300 data points / spectrum
wavelength: PC [integer] 1 2 ... 300
data: (875 rows x 4 columns)
  1. y: y/(mu * m) [numeric] -4.77 -4.77 ... 19.23
  2. x: x/(mu * m) [numeric] -11.55 -10.55 ... 22.45
  3. clusters: clusters [factor] matrix matrix ... lacuna + NA
  4. spc: score / a.u. [AsIs matrix x 300] 1.1892 1.3720 ... 3.7274e-17
```

The loadings can be similarly re-imported:

```
> loadings <- decomposition (chondro, t(pca$rotation), scores = FALSE,
+                             label.spc = "loading I / a.u.")
> loadings
```

```
hyperSpec object
  300 spectra
   1 data columns
  300 data points / spectrum
wavelength: Delta * tilde(nu)/cm^-1 [numeric] 602 606 ... 1798
data: (300 rows x 1 columns)
  1. spc: loading I / a.u. [AsIs matrix x 300] 0.016709 0.007826 ... -0.092112
```

There is, however, one important difference. The loadings are thought of as values computed from all spectra together. Thus no meaningful extra data can be assigned for the loadings object (at least not if the column consists of different values). Therefore, the loadings object lost all extra data (see above).

*retain.columns* triggers whether columns that contain different values should be dropped. If it is set to `TRUE`, the columns are retained, but contain NAs:

```
> loadings <- decomposition (chondro, t(pca$rotation), scores = FALSE,
+                             retain.columns = TRUE, label.spc = "loading I / a.u.")
> loadings[1]$..
```

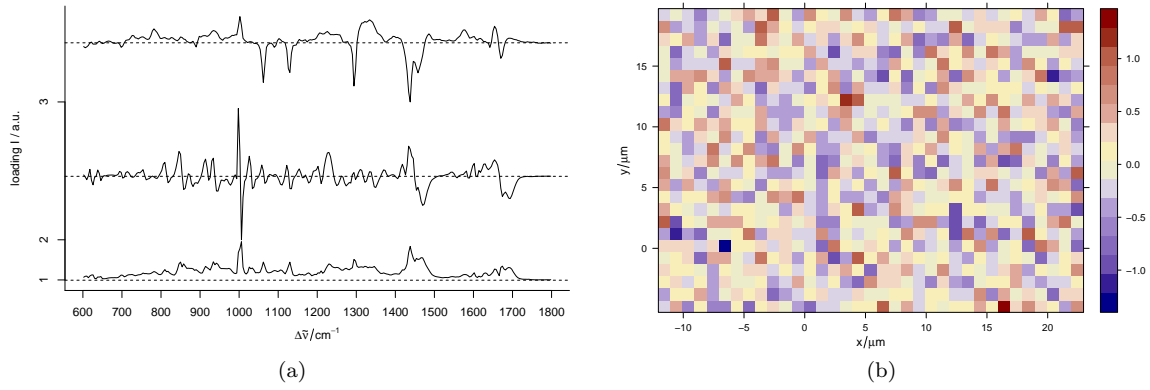


Figure 6: (a) The first three loadings: `plot (loadings [1 : 3], stacked = TRUE)`. (b) The second score map: `plotmap (scores [, , 2])`.

```

      y x clusters
1 NA NA      <NA>

```

If an extra data column does contain only one unique value, it is retained anyways:

```

> chondro$measurement <- 1
> loadings <- decomposition (chondro, t(pca$rotation), scores = FALSE,
+                             label.spc = "loading I / a.u.")
> loadings[1]$..

      measurement
1                1

```

### 13.1.1. PCA as Noise Filter

Principal component analysis is sometimes used as a noise filtering technique. The idea is that the relevant differences are captured in the first components while the higher components contain noise only. Thus the spectra are reconstructed using only the first  $p$  components.

This reconstruction is in fact a matrix multiplication:

$$spectra^{(nrow \times nwl)} = scores^{(nrow \times p)} loadings^{(p \times nwl)}$$

Note that this corresponds to a model based on the Beer-Lambert law:

$$A_n(\lambda) = c_{n,i} \epsilon(i, \lambda) + error$$

The matrix formulation puts the  $n$  spectra into the rows of  $A$  and  $c$ , while the  $i$  pure components appear in the columns of  $c$  and rows of the absorbance coefficients  $\epsilon$ .

For an ideal data set (constituents varying independently, sufficient signal to noise ratio) one would expect the principal component analysis to extract something like the concentrations and pure component spectra.

If we decide that only the first 10 components actually carry spectroscopic information, we can reconstruct spectra with better signal to noise ratio:

%%

```

> smoothed <- scores [, , 1:10] %*% loadings [1:10]


```

Keep in mind, though, that we cannot be sure how much *useful* information was discarded with the higher components. This kind of noise reduction may influence further modeling of the data. Mathematically speaking, the rank of the new  $875 \times 300$  spectra matrix is only 10.

### 13.2. Data Analysis using long-format data.frame e.g. plotting with ggplot2

Some functions need the data being an *unstacked* or *long-format data.frame*. `as.long.df` is the appropriate conversion function. `as.long.df`

```
> library (ggplot2)
> p <- ggplot (as.long.df (chondro [1]), aes (x = .wavelength, y = spc)) + geom_line ()
```

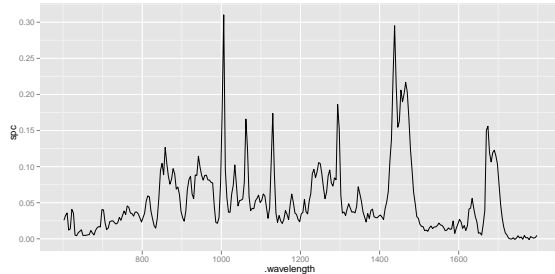

### 13.3. Data Analysis Methods using a matrix e.g. Hierarchical Cluster Analysis

`[[]]`

```
> dist <- pearson.dist (chondro [[]])
> dendrogram <- hclust (dist, method = "ward")

> plot (dendrogram)
```

In order to plot a cluster map, the cluster membership needs to be calculated from the dendrogram. First, cut the dendrogram so that three clusters result:

```
> chondro$clusters <- as.factor (cutree (dendrogram, k = 3))
```

As the cluster membership was stored as factor, the levels can be meaningful names, which are displayed in the color legend.

```
> levels (chondro$clusters) <- c ("matrix", "lacuna", "cell")
```

Then the result may be plotted (figure 7b):

### 13.4. Calculating group-wise Sum Characteristics e.g. Cluster Mean Spectra

`aggregate` applies the function given in *FUN* to each of the groups of spectra specified in *by*.

`aggregate`

So we may plot the cluster mean spectra:

```
> means <- aggregate (chondro, by = chondro$clusters, mean_pm_sd)
> plot (means, col = cluster.cols, stacked = ".aggregate", fill = ".aggregate")
```

### 13.5. Splitting an Object, and Binding a List of hyperSpec Objects

`split`

A *hyperSpec* object may also be split into a list of *hyperSpec* objects:

```
> clusters <- split (chondro, chondro$clusters)
> clusters
```

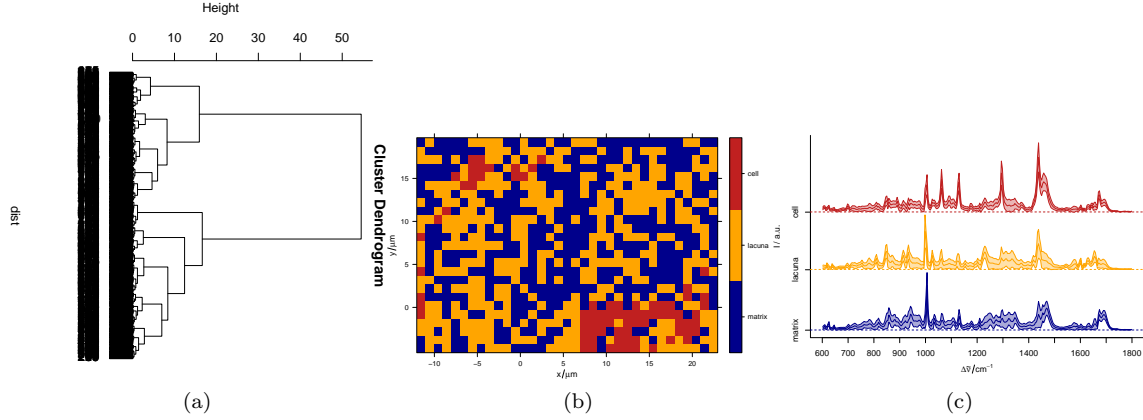


Figure 7: The results of the cluster analysis: (a) the dendrogram (b) the map of the 3 clusters (c) the mean spectra.

```
$matrix
hyperSpec object
  404 spectra
  5 data columns
  300 data points / spectrum
wavelength: Delta * tilde(nu)/cm^-1 [numeric] 602 606 ... 1798
data: (404 rows x 5 columns)
  1. y: y/(mu * m) [numeric] -4.77 -4.77 ... 19.23
  2. x: x/(mu * m) [numeric] -11.55 -9.55 ... 22.45
  3. clusters: clusters [factor] matrix matrix ... matrix
  4. spc: I / a.u. [matrix300] 0.025669 0.040389 ... 0.0034843
  5. measurement: measurement [numeric] 1 1 ... 1

$lacuna
hyperSpec object
  392 spectra
  5 data columns
  300 data points / spectrum
wavelength: Delta * tilde(nu)/cm^-1 [numeric] 602 606 ... 1798
data: (392 rows x 5 columns)
  1. y: y/(mu * m) [numeric] -4.77 -4.77 ... 19.23
  2. x: x/(mu * m) [numeric] -10.55 -8.55 ... 20.45
  3. clusters: clusters [factor] lacuna lacuna ... lacuna
  4. spc: I / a.u. [matrix300] 0.029119 0.044025 ... 0.0009416
  5. measurement: measurement [numeric] 1 1 ... 1

$cell
hyperSpec object
  79 spectra
  5 data columns
  300 data points / spectrum
wavelength: Delta * tilde(nu)/cm^-1 [numeric] 602 606 ... 1798
data: (79 rows x 5 columns)
  1. y: y/(mu * m) [numeric] -4.77 -4.77 ... 17.23
  2. x: x/(mu * m) [numeric] 7.45 8.45 ... 2.45
  3. clusters: clusters [factor] cell cell ... cell
  4. spc: I / a.u. [matrix300] 0.045152 0.049662 ... -8.4008e-06
  5. measurement: measurement [numeric] 1 1 ... 1
```

Splitting can be reversed by `rbind` (see section 9.1, page 8). Another, similar way to combine a number of `hyperSpec` objects with different wavelength axes or extra data columns is `collapse` (see section 9.2, page 9).



## 14. Speed Considerations

While most of *hyperSpec*'s functions work at a decent speed for interactive sessions (of course depending on the size of the object), iterated (repeated) calculations as for bootstrapping or iterated cross validation may ask for special speed considerations.

In that case, the calculations may be sped up considerably if the required parts of the *hyperSpec* object are extracted into a *data.frame* or matrices beforehand. A related tip is that many model fitting functions in R are much faster if the formula interface is avoided and the appropriate *data.frames* or matrices are handed over directly.

Another possibility to speed up is to switch of the automatic logging of how the objects are transformed. Logging involves appending rows to the *data.frame* in slot `@log`. While the absolute amount of time needed to add a logbook entry is small, it may be executed very often (e.g. during each call of `[]`). However, the first strategy will usually allow for far larger gains in speed.

As an example, let's again consider the code for shifting the spectra:

```
> system.time ({
+   for (i in seq_len (nrow (chondro)))
+     chondro [[i]] <- interpolate (chondro [[i]], shifts [i], wl = wl (chondro))
+ })

   user system elapsed
43.490   1.270  44.913

> rm (chondro)

> hy.setOptions (log = FALSE)
> system.time ({
+   for (i in seq_len (nrow (chondro)))
+     chondro [[i]] <- interpolate (chondro [[i]], shifts [i], wl = wl (chondro))
+ })

   user system elapsed
28.290   1.330  29.701

> hy.setOptions (log = TRUE)
> rm (chondro)

> system.time ({
+   tmp <- chondro [[]]
+   for (i in seq_len (nrow (chondro)))
+     tmp [i, ] <- interpolate (tmp [i, ], shifts [i], wl = wl (chondro))
+   chondro [[]] <- tmp
+ })

   user system elapsed
 0.920   0.000   0.924

> rm (chondro)
```

## References

- [1] Ron Wehrens and Bjørn-Helge Mevik. *pls: Partial Least Squares Regression (PLSR) and Principal Component Regression (PCR)*, 2007. URL <http://mevik.net/work/software/pls.html>. R package version 2.1-0.

## A. Overview of the functions provided by *hyperSpec*

Function	Explanation
<i>Access parts of the object</i>	
[	Select / extract / delete spectra, wavelength ranges or extra data
[<-	Set parts of spectra or extra data
[[	Select / extract / delete spectra, wavelength ranges or extra data, get the result as matrix or data.frame
[[<-	Set parts of spectra matrix
\$	extract a data column (including \$spc)
\$<-	replace a data column (including \$spc)
i2wl	convert spectra matrix column indices to wavelengths
isample	get a random sample of the spectra as index vector
labels	get column labels
labels<-	set column labels
logbook	logging the data treatment
logentry	make a logbook entry
rownames<-	
sample	generate random sample of the spectra
seq.hyperSpec	sequence along the spectra, either as <i>hyperSpec</i> object or index vector
wl	extract the wavelengths
wl<-	replace the wavelengths
wl2i	convert wavelengths to spectra matrix column indices
<i>Basic information</i>	
colnames	
colnames<-	
dim	
dimnames	
length	
ncol	number of data columns (extra data plus spectra matrix)
nrow	number of spectra
nwl	number of data points per spectrum
print	summary information
rownames	
show	
summary	summary information including the log
chk.hy	checks whether the object is a hyperSpec object
<i>Combine/split</i>	
bind	commom interface for rbind and cbind

Function	Explanation
<code>cbind2</code>	bind two <i>hyperSpec</i> objects by column
<code>cbind.hyperSpec</code>	
<code>collapse</code>	combine objects by adding columns if necessary. See <code>plyr::rbind.fill</code> .
<code>rbind2</code>	bind two <i>hyperSpec</i> objects by row, i. e. add wavelength ranges or extra data
<code>rbind.hyperSpec</code>	bind objects by row, i. e. add wavelength ranges or extra data
<code>split</code>	
<code>merge</code>	combines spectral ranges. works if spectra are in only one of the data sets
<i>Comparison</i>	
<code>all.equal</code>	
<code>Compare</code>	<code>&gt; &lt; == &gt;= &lt;=</code> return a logical matrix
<code>is.na</code>	
<i>Create and initialize an object</i>	
<code>initialize</code>	
<i>File import/export</i>	
<code>read.ENVI</code>	import ENVI file
<code>read.ENVI.Nicolet</code>	import ENVI files written by Nicolet spectrometers
<code>read.spc</code>	import .spc file
<code>read.spc.KaiserMap</code>	import a Raman map saved by Kaiser Optical Systems' Hologram software as multiple .spc files
<code>read.txt.long</code>	import long-type ASCII file
<code>read.txt.wide</code>	import wide-type ASCII file
<code>R.matlab::readMat</code>	import matlab file
<code>R.matlab::writeMat</code>	export as matlab file
<code>scan.txt.Renishaw</code>	import ASCII files produced by Renishaw (InVia) spectrometers
<code>write.txt.long</code>	export as long-type ASCII file
<code>write.txt.wide</code>	export as wide-type ASCII file
<code>scan.zip.Renishaw</code>	directly read zip packed ASCII files produced by Renishaw spectrometers
<i>Maths</i>	
<code>%*%</code>	matrix multiplication
<code>Arith</code>	
<code>log</code>	
<code>Math</code>	mathematical functions. See <code>(help ("Math extquotedbl ))</code>
<code>Math2</code>	rounding
<code>Summary</code>	summary measures such as <code>min</code> , <code>max</code> , etc.

Function	Explanation
<i>Plotting</i>	
<code>levelplot</code>	
<code>map.identify</code>	identify spectra in map plot
<code>matlab.dark.palette</code>	darker version of <code>matlab.palette</code>
<code>matlab.palette</code>	palette resembling Matlab's jet colors
<code>plot</code>	main switchyard for plotting
<code>plotc</code>	intensity over one other dimension: calibration plots, time series, depth series, etc.
<code>plotmap</code>	false-colour intensity over two other dimensions: spectral images, maps, etc. (rectangular tessellation)
<code>plotspc</code>	spectra plots: intensity over wavelength
<code>plotvoronoi</code>	false-colour intensity over two other dimensions: spectral images, maps, etc. (Voronoi tessellation)
<code>spc.identify</code>	identify spectra and wavelengths in spectra plot
<code>stacked.offsets</code>	calculate intensity axis offsets for stacked spectral plots
<code>trellis.factor.key</code>	modify list of <code>levelplot</code> arguments according to factor levels
<i>Spectra-specific transformations</i>	
<code>orderwl</code>	sort columns of spectra matrix according to the wavelengths
<code>spc.bin</code>	spectral binning
<code>spc.fit.poly</code>	least squares fit of a polynomial
<code>spc.fit.poly.below</code>	least squares fit of a polynomial with automatic support point determination
<code>spc.loess</code>	<code>loess</code> smoothing interpolation
<i>Type conversion</i>	
<code>as.character</code>	
<code>as.data.frame</code>	
<code>as.long.df</code>	convert to a long-format data.frame.
<code>as.matrix</code>	
<code>as.wide.df</code>	convert to a wide-format data.frame with each wavelength one column
<code>decomposition</code>	re-import results of spectral matrix decomposition (or the like) into <i>hyperSpec</i> object
<i>Utility functions</i>	
<code>array2df</code>	convert array into a matrix or data.frame
<code>array2vec</code>	convert array indices ( $n$ element vector) into vector indices
<code>mean_pm_sd</code>	mean $\pm$ one standard deviation of a vector
<code>mean_sd</code>	mean and standard deviation of a vector
<code>pearson.dist</code>	distance measure based on Pearson's $R^2$
<code>rbind.fill</code>	transitional patch of <code>plyr::rbind.fill</code> working with matrices

Function	Explanation
<code>rbind.fill.matrix</code>	transitional until <code>plyr::rbind.fill.matrix</code> is out
<code>vec2array</code>	convert vector (one element) index into an array into an $n$ element array index
<code>wc</code>	word count using <code>wc</code> if available on the system
<i>Vectorization</i>	
<code>aggregate</code>	
<code>apply</code>	
<code>sweep</code>	

## Session Info

R version 2.12.1 (2010-12-16)

Platform: x86\_64-pc-linux-gnu (64-bit)

locale:

```
[1] LC_CTYPE=de_DE.utf8      LC_NUMERIC=C             LC_TIME=de_DE.utf8
[4] LC_COLLATE=de_DE.utf8    LC_MONETARY=C            LC_MESSAGES=de_DE.utf8
[7] LC_PAPER=de_DE.utf8      LC_NAME=C                LC_ADDRESS=C
[10] LC_TELEPHONE=C           LC_MEASUREMENT=de_DE.utf8 LC_IDENTIFICATION=C
```

attached base packages:

```
[1] grid      stats      graphics  grDevices utils      datasets  methods   base
```

other attached packages:

```
[1] ggplot2_0.8.9      proto_0.3-8      reshape_0.8.4      plyr_1.4
[5] plotrix_3.0-8      hyperSpec_0.96-20110208 lattice_0.19-17
```